

**1300 LPM
TRUE AIR
LOSS
BLOWER
SYSTEM!**

NEW!



Available
1 Pack

**Future Air True Low Air Loss
Mattress System**

14000 79" (L) x 35.5" (W) x 10" (H), 1/ea

HCPCS code: E0277

FEATURES:

- True low air loss, immersion and envelopment therapies to optimally manage skin maceration, pressure, shear and friction.
- Zoned therapy customizes bladder inflation to 3 zones to maximize immersion and redistribute pressure around bony prominences.
- Automatic fowler boost when head of bed is elevated.
- Segmented side bolsters allow for safe and easy bed frame articulation.
- 10" side bolster provides safety support.
- Auto 3-minute panel lockout prevents accidental setting change.
- Max Firm mode helps facilitate nursing care.
- Quick connector auto-detects ON/OFF and facilitates CPR deflation.
- 440 lb. weight capacity.
- Limited 12 Month Warranty on pump and mattress.



Heavy Duty
Versatile Hanger




Auto Boost
Sensor Pad

**3 DIFFERENTIAL
PRESSURE ZONES
FOR HEAD, BODY
AND FOOT SECTIONS
TO PROVIDE
CUSTOMIZED
PRESSURE
REDISTRIBUTION.**



Quick Disconnect

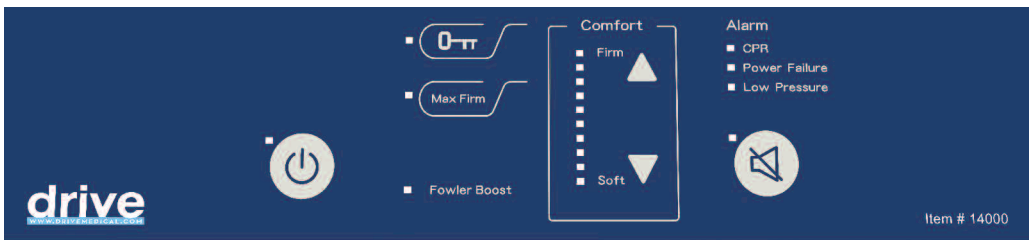

SPECIFICATIONS

Electrical	AC 100-120V 50/60 Hz, 5.0A
Pressure	Adjustable 10 level pressure range
Air Output (lpm)	Maximum airflow 1,300 LPM
PUMP Dimensions (L x W x H)	6.2" x 15" x 10.5"
PUMP Weight	12.1 lb.
Mattress Model	8" Mattress with 10" side bolster 2" Foam pocket is provided
Cover Mattress	4-Way stretch, vapor permeable, low shear cover
Mattress Weight	15.4 lb.
Noise Level	46 dBA (At comfort level 5, measured 1 meter from pump)
Weight Capacity	440 lb.

pressure prevention

To order call toll free: 877.224.0946

drive