

## Med-Aire Plus Bariatric Alternating Pressure Mattress Replacement System

14054 80" (L) x 54" (W) x 10" (H), Pump and Mattress System, 1/ea

**HCPCS code: E0277**

### FEATURES:

- This bariatric system provides alternating pressure to optimize pressure redistribution for patients with up to Stage IV pressure ulcers.
- Designed to prevent, treat and manage pressure ulcers in the home, long term care, or acute care setting.
- **Cell-on-cell design** of 10" deep air cells prevents "bottoming out" and also provides up to 24 hours of power outage protection.
- **Digital pump** allows for customization to individual patients' comfort and safety needs, and will auto return to alternating pressure mode if system is left on static mode for >2 hours.
- **4-way stretch cover** is low shear, fluid-resistant, vapor permeable, quilted and zippered.
- Each one of the 20, **10" deep air bladders** are easily removed and replaced.
- Pillow function maintains the 3 air cells at the head of the bed in static mode for patient comfort.
- 3 diagnostic alarms tell you whether a problem is with the mattress or the pump for more efficient servicing.
- 1000 lb. weight capacity.
- CPR valve allows for rapid deflation.



14054



14427 (Electric Pump)

Available  
**1 Pack**

**54" WIDE!**

BARIATRIC

**1000 LB**

WEIGHT CAPACITY

### PARTS:

14054M	80" (L) x 54" (W) x 10" (H), Mattress Only, 1/ea
14030XP	Pump Only for 14054, 1/ea
14054-CELL	Replacement Bladder for 14054, 1/ea

### Mattress:

- Dimensions: 80" (L) x 54" (W) x 10" (H)
- Base Material: zippered vinyl
- **Maximum Support Capacity: 1000 lbs.**
- Limited 18 Month Warranty

### SPECIFICATIONS:


#### Pump:

- Power: 120 VAC, 60 Hz
- Pressure Range: 20 mmHg to 60 mmHg (inner mattress pressure)
- Auto firm (60 mmHg, Auto return in 30 min.)
- Seat function (Pressure preset + 10 mmHg)
- Static selectable with auto-return to alternating mode.
- Adjustable cycle time 10, 15, 20, 25 minutes.
- Function lock (Auto lock in 30 seconds)
- Low pressure alarm
- System (alternating) alarm
- Power failure alarm
- Alarm mute
- Memory function
- **Air Flow: 12-LPM**
- Power Cord Length: 15 Ft.
- Pump Dimensions: 4.5" (L) x 14.2" (W) x 9.8" (H)
- UL: 60601-1, 1st Edition-2003
- Limited 18 Month Warranty

**NEW  
DIGITAL  
PUMP  
DESIGN!**



14030XP  
**PUMP ONLY  
NEW & IMPROVED!**

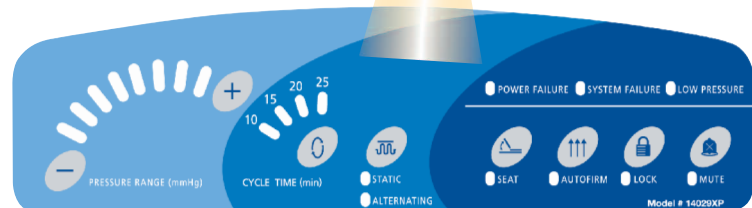

\*Important Notice: Item numbers that display an asterisk (\*) reflect an item that will be discontinued or package change.