



Find Out More on the Internet  
[www.harmaraccess.com](http://www.harmaraccess.com)




Still want to know more about Harmar Access and our products? Visit us on the web. You'll find more photos and detailed information.  
[www.HarmarAccess.com](http://www.HarmarAccess.com)

**More Harmar Products**

Harmar Access produces Bath Lifts, Porch Lifts and Stair Lifts that enhance peoples mobility independence and quality of life. We stand committed to bring people the highest quality Access and Mobility products.

**Stair Lifts**



Harmar Access Stair Lifts transform unfriendly multi-story houses into easily accessible homes. Models available for every budget make staying in your home a safe and simple reality. Harmar Access Stair Lifts are the perfect solution to reliably and effortlessly traverse your staircase in complete comfort.

**Porch Lifts**



Harmar Access Porch Lifts are safe, smooth and economical solutions to the barriers stairs can create outside a home. Available in 2, 4 and 6-foot vertical lifting heights, our residential lifts are designed to be simple to install and operate and are virtually maintenance free for use indoors and out.

**Harmar Mobility**



Harmar Mobility produces lift and ramp products that enhance people's mobility, independence and quality of life. Products include Inside and Outside Lift Systems, Hybrid Platform Lift Systems and Portable and Residential Ramp Systems. [www.harmar.com](http://www.harmar.com)



[www.harmaraccess.com](http://www.harmaraccess.com)



AquaJoy  
 made by :

[www.harmaraccess.com](http://www.harmaraccess.com)

**The UK's  
 Number  
 One Bathlift**

**H**armar Access is the exclusive North American distributor for the leading European bathtub lift, the AquaJoy. With no installation required, the AquaJoy Premier PLUS Bath Lift makes getting into and out of the bathtub safe, simple and easy. Its universal design fits every standard bathtub and features automatic reclining as well as a low battery safety switch. As with all Harmar products, the AquaJoy Premier PLUS continues the tradition of high quality products designed to be easy to use and simple to install.

Harmar Access builds upon our commitment to mobility by breaking down barriers in the home and making every room in the house accessible. Harmar Access has a full line of Porch Lifts, Stair Lifts and Bathroom Lifts designed to safely increase your mobility and home accessibility.

**Chad C. Williams**  
 President, Harmar Access



## AquaJoy Premier PLUS Bathlift

The AquaJoy Premier PLUS is a two piece, lightweight, battery operated bathlift which fits into any style of bath – even corner baths. It can be electronically reclined or stopped at the desired sitting position, allowing the individual to choose the most comfortable bathing position for them.

The bathlift maximizes legroom in all baths and with the easy release suction feet and lifting bar as standard, the Premier PLUS can be taken in and out of the bath with great ease. It is the easiest bathlift to put together and take apart, ensuring ease of use and transportation for all involved. The extra wide, comfortable seat allows the user to be stable and safe at all times.

The AquaJoy Premier PLUS can be used with or without cover mats. The standard cover mats are easily removed and can be washed in any normal washing machine up to 140 degrees.

The fully waterproof, floating, lightweight operating hand control operates easily. With convex and concave buttons as well as an indicator light, the controller is ideal for all types of users. The battery safety mechanism will not allow the bathlift to lower unless there is sufficient charge to raise the bathlift back to the top of the bath. The AquaJoy lifting capacity ensures the bathlift will also be available to clients with weight lifting capabilities up to 375 pounds.

The AquaJoy range of accessories including gel Swivel Slide, Headrest with Wings, Chest Harness with Pommel, and Lap Harness with Pommel can be added onto any Premier PLUS bathlift, allowing versatility for clients and better options for Loan Equipment Services. The AquaJoy Premier PLUS is produced with NANO technology to reduce the risk of cross contamination and infection. Sideflap / Handle diverters are packed free of charge in every box and additional suction cups can be added for baths with awkward ripples or dimples.

**Easy to set up!**




### AquaJoy Key Features

- \*Comfortable, washable padding
- \*Automatically reclines up to 40 degrees (user selected)
- \*Battery operated with low-battery safety switch
- \*Suction cup feet for stability
- \*Waterproof, floating hand controller with large buttons
- \*Simple to disassemble for storage and transportation
- \*375lb. weight capacity
- \*Easy to use battery charger included
- \*Many accessories available

### AquaJoy Specifications

- \*Length of seat: 22.8 inches
- \*Width of seat with extended seat side flaps: 27 inches
- \*Width of seat with seat side flaps folded in: 15.3 inches
- \*Height of seat in lowest position: 3.1 inches
- \*Height of seat in its highest position: 18.1 inches
- \*Depth of seat measured from Backrest to front seat edge: 23.8 inches
- \*Height of backrest: 23.8 inches
- \*Width of backrest: 13.4 inches
- \*Carrying capacity: 375 lb.
- \*Total Weight without battery: 27.5 lb.
- \*Weight of seat: 16.3 lb.
- \*Weight of backrest: 11.2 lb.
- \*Weight of battery hand controller: 1.4 lb.
- \*Duration of stroke (176 lb. load): Up approx. 50 sec. Down approx. 35 sec.
- \*Charging time with flat battery: approx. 4hours

### AquaJoy Bath Lift Accessories

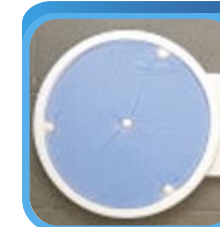

#### BL105 Swivel Slide

Eases entry and exit from your bath. Our most popular accessory



#### BL120 Headrest with Wings

Provides extra support and security for the head and upper body.



#### BL115 Chest Harness with Pommel

Ensures occupant safety by preventing downward slippage.



#### BL110 Lap Harness with Pommel

Ensures occupant safety by preventing any slippage for patients with compromised upper body strength.