

MEDTRADE EDITION

P O W E R   C H A I R



CHERRY RED



PARADISE BLUE



SOLAR FLARE ORANGE  
*Limited Edition*



ENVY GREEN  
*Limited Edition*



SATIN SILVER  
*Limited Edition*



SUNBURST YELLOW  
*Limited Edition*

Seat Color



Black Vinyl




*Available in 6 New  
Vibrant Colors with  
a Larger Foot Plate and  
New Shroud Design!*

**Disassembles  
Quickly &  
Easily!**

LiteRider<sup>®</sup> GOLDEN  
ENVY  
2016

The all new LiteRider Envy Power Wheelchair is both stylish and durable with a tubular design that is unique to the Golden LiteRider family! The LiteRider Envy, featuring 2-22 amp hour batteries for an operating range of up to 15.5 miles, is Golden's Retail Cash Power Wheelchair! The LiteRider Envy has a larger footplate and an impressive turning radius of 28.5", and like the other LiteRiders in the family, it disassembles quickly and easily for transport. The heaviest piece is only 37 lbs., making storage convenient in the trunk of any vehicle. When you are looking for a compact power wheel chair that is fun to drive and easy to transport – look no further than the Golden LiteRider Envy!



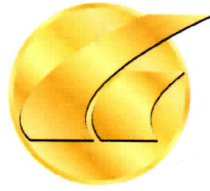

# INTRODUCING THE NEW

## LiteRider<sup>®</sup> GOLDEN *ENVY* 2018

America's favorite full-time portable power chair now available in 6 New Vibrant Colors with a Larger Foot Plate and New Shroud Design!

**Ships Now:** Cherry Red & Paradise Blue

**Ships in Early 2018:**  
Satin Silver, Solar Flare Orange,  
Envy Green & Sunburst Yellow



LIMITED EDITION  
**SHOWN IN ENVY GREEN**

**CHERRY RED**



LIMITED EDITION  
**SUNBURST YELLOW**



LIMITED EDITION  
**SOLAR FLARE ORANGE**



**PARADISE BLUE**



LIMITED EDITION  
**SATIN SILVER**





### Model Number

Medicare Code  
Weight Capacity  
Drive Wheels  
Maximum Speed  
Operating Range<sup>1</sup>  
Ground Clearance  
Turning Radius  
Type Batteries

### Colors

Freewheel Mode  
Electronic Speed Control  
Electro-Mechanical Brakes  
Charger  
Controller

### Dimensions:

Length  
Width  
Height (ground to top of back of seat)  
Height (ground to top of back of headrest)  
Ground to Top of Seat  
Deck to Top of Seat  
Number of Seat Height Positions  
Size of Increments

### Weight of Unit: (assembled)

Front Half  
Rear Half  
Seat w/Arms  
Batteries

### Tires:

Drive Wheels  
Casters  
Anti-Tip Wheels

### Standard Seat:

Back Height - no headrest  
Back Height with headrest  
Width x Depth  
Color

### Optional Seats:

Back Height - no headrest  
Back Height with headrest  
Width x Depth  
Color

### Seat Options:

Pneumatic Rise Seat Option  
Electric Rise Seat Option

### Footrest:

Height adjustable  
Angle adjustable

### Warranty:

Frame  
Drive Train  
Electronics<sup>2</sup>  
Batteries

### GP162

Not Coded  
300 lbs.  
Rear  
3.5 mph  
15.5 miles  
2.5"  
28.5"  
2-22 AH  
Cherry Red, Paradise Blue,  
Solar Flare Orange,  
Envy Green, Satin Silver,  
Sunburst Yellow  
Yes  
Yes  
Yes  
Off Board  
PG VR-2  
50 Amp Remote

36"  
22.5"  
33" - 34"  
N/A  
20.5" - 21.5"  
N/A

123 lbs.  
37 lbs.  
36 lbs.  
23 lbs.  
27 lbs/pack

Flat Free  
9"  
6"  
3"

Stadium Style  
16"  
N/A  
20" x 18"  
Black Vinyl  
Stadium Style  
16"  
N/A  
20" x 18"  
Black Vinyl

N/A  
N/A  
Yes  
Yes  
Yes

Five Years  
13 Month Warranty  
13 Month Warranty  
12 Month Warranty

## Disassembles easily for transport!




1. Remove seat



2. Pull out battery pack



3. Grasp the frame and locking handle; lift the handle to unlock, and lift the frame up and off the drive train.

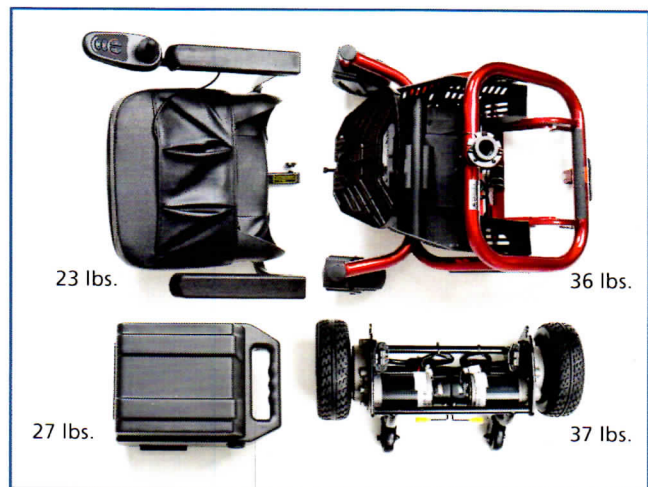

Perfect for transport!

## Available Accessories:

- Cane Holder Tube
- Oxygen Tank Holder
- Cup Holder
- Quad Cane Holder
- Walker Holder
- Crutch Holder
- Safety Flag
- Power Chair Cover
- Pack N' Go
- Seat Belt (Optional)

Some accessories require a mounting bracket or mounting clips.


<sup>1</sup> Battery range tested at 200 lbs. and will vary due to rider weight, drive surface and terrain.

<sup>2</sup> Electronics warranty excludes batteries.



# BUZZAROUND EX CHARGING DOCK INSTRUCTIONS



 **WARNING INDOOR USE ONLY!**

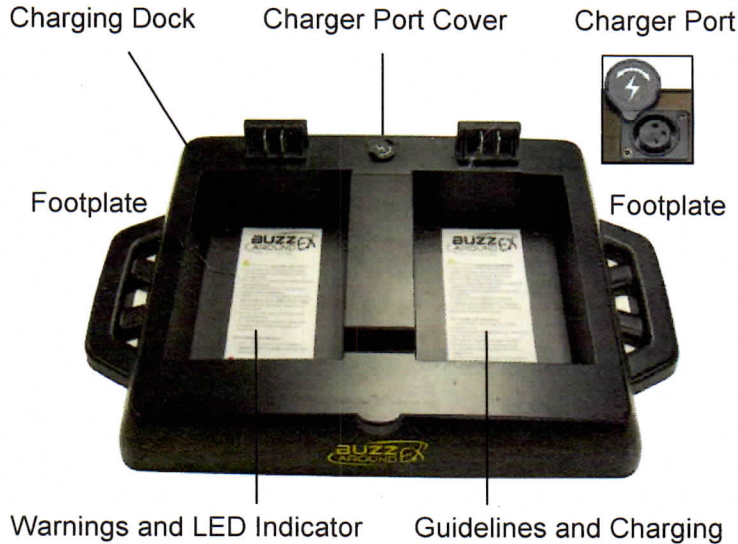

Battery Packs Installed - Ready to Charge


**NOTICE** Before attempting to charge your two battery packs , please read and follow all warnings, guidelines, and instructions affixed to the product.

## Fuse Replacement

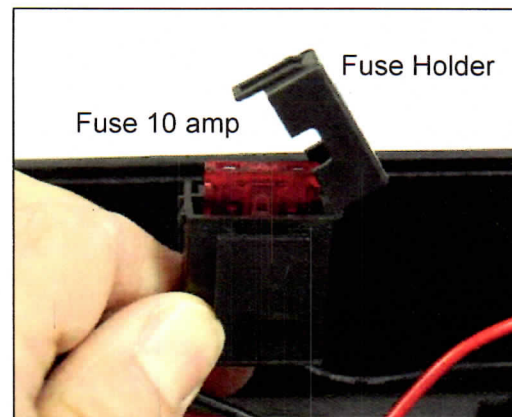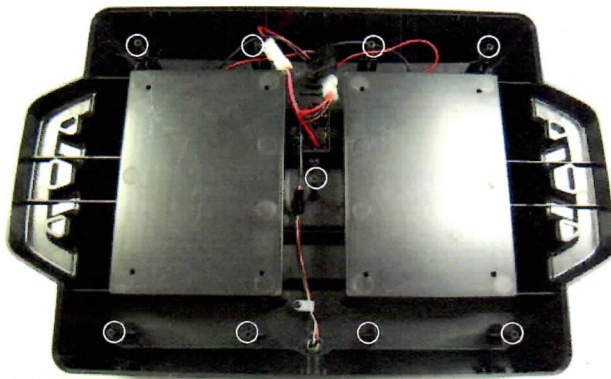


Fuse Holder  
Fuse 10 amp


1. Remove the nine screws on the bottom side of the charging dock to remove the cover.
2. Open the fuse holder.
3. Remove the (blade style) 10 amp fuse.
4. Reverse steps 1-3 to install the new fuse.

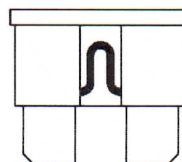

Good Fuse

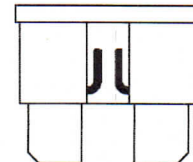

Bad Fuse

 **WARNING Fuse must be replaced with the exact type and rating.**