


AL580 & AL580XL

Next Generation Power Wheelchair Lifts



Lift folds automatically when not in use

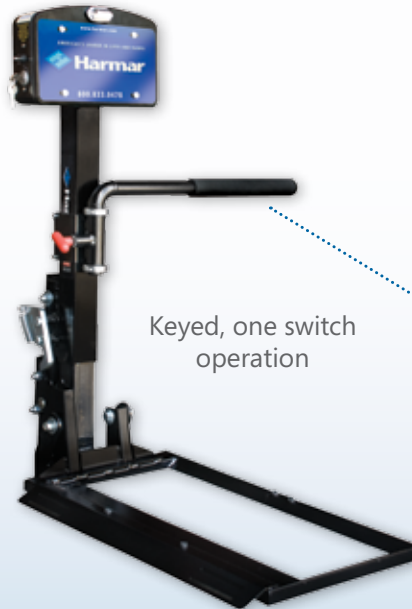

Keyed, one switch operation



Adjustable padded hold down arm provides a secure, custom fit for every chair without the use of straps or tie-downs

Lightweight, unobtrusive platform designed for power wheelchairs

FEATURES

- Padded hold-down arm conveniently secures power wheelchair
- No power wheelchair modifications required
- Lift folds automatically when unoccupied
- No straps or tie-downs required
- Manual crank backup & license plate mount included
- Swing-away option available

SPECIFICATIONS

Lifting Capacity:	AL580 - 350 lbs AL580XL - 350 lbs
Platform Size:	AL580 - 16.5" x 27" AL580XL - 17.5 x 27"
AL580 Designed For:	MWD Power Chair with 10" Wheels
AL580XL Designed For:	MWD Power Chair with 14" Wheels
Installed Weight:	83 lbs
Hitch Height (range):	5"-24" <small>(With appropriate hitch adapter installed)</small>
Hitch Class:	AL580 - Class II or III AL580XL - Class III
Warranty:	3-Year Transferable

DESCRIPTION

The lightest solutions on the market providing more options! Compact, durable and American-made, the AL580 and AL580XL outside vehicle platform lifts safely and securely transport power wheelchairs. With a tidy and space efficient rear profile, the AL580 and AL580XL have a highly desirable and understated presence on the vehicle exterior.

The AL580 and AL580XL are easy to operate and are built to last for years of trouble-free use. Designed to meet the needs of consumers who need a quick, automatic securing option without making modifications to their power wheelchair, the stylish and compact detailing of this lift make it a must have for power wheelchair owners.

AL580 & AL580XL Next Generation Power Wheelchair Lifts



Only 83 lbs Installed!

AL580 / AL580XL



115 lbs Installed!

The 'Other' Guys

In many cases these lifts are the only lift solution on the market that will work. While other companies are raising prices and cheapening materials, Harmar is working to bring newer and more innovative products to market.

Easy as 1, 2, 3!



1

Drive chair onto platform



2

Raise armrest closest to vehicle and rotate swing arm out over top of power wheelchair



3

Raise lift and go!