



LIGHTWEIGHT PORTABLE SCOOTER

Body Colors



Red



Blue



Sahara

Improved!
Delta Tiller
and larger seat
now standard!

Ride

BUZZaround[®]
Lite

3 wheel:
 GB-106 & GB-106XR

The popular Buzzaround Lite scooter has been improved! This compact scooter now comes standard with a Delta tiller with wraparound handle and a larger, more comfortable seat. The Buzzaround Lite still features our easy "Pop! Pull! Push!" wireless disassembly for transport. Change the look of your scooter with our three sets of snap-on shroud panels in Red, Blue and Sahara, all while enjoying its tight indoor maneuverability and smooth outdoor use.



Rise.

Roll.

Rest.

Ride.



GOLDEN
 technologies

YOUR LIFE IN MOTION

3-wheel Specifications



Model	GB-106	GB-106XR
Length	36.5"	36.5"
Width	21.25"	21.25"
Tires	Foam-filled	Foam-filled
Front tire	8"	8"
Rear tire	8"	8"
Max speed	4 mph	4 mph
Ground clearance	2"	2"
Battery range using 2-12Ah batteries At maximum weight (250 lbs.)	8 miles*	
Battery range using 2-15Ah batteries At maximum weight (250 lbs.)		10 miles*
Turning radius	31"	31"
3 point turn	39"	39"
90 degree corridor turn	28"	28"
Seat weight with arms	21 lbs.	21 lbs.
Battery pack weight		
2 x 12 Amp hour batteries	21 lbs.	
2 x 15 Amp hour batteries		23 lbs.
Front half weight	32 lbs.	32 lbs.
Rear half weight	24 lbs.	24 lbs.
Freewheel mode	Yes	Yes
Electronic speed control	Yes	Yes
Brakes:	Electro-mechanical	Electro-mechanical
Charger	Off board	Off board
Rear-wheel drive	sealed transaxle, motor 24 VDC	sealed transaxle, motor 24 VDC
Controller	45 amp/programmable	45 amp/programmable
Warranty – frame	Limited lifetime	Limited lifetime
Warranty – electronics**	2 years	2 years
Warranty – drive train	2 years	2 years

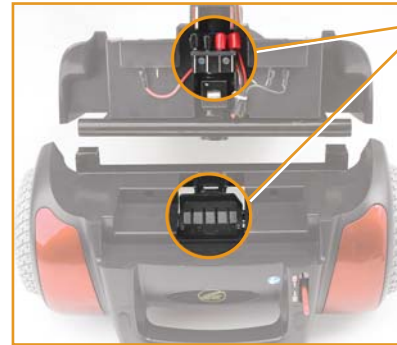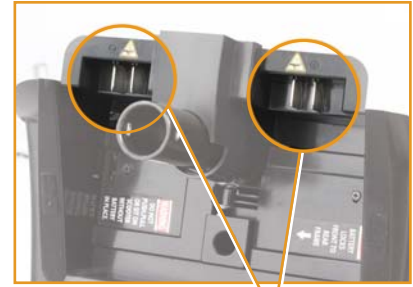

Exclusive interlocking teeth provide the tightest, most dependable battery connection for lasting power.




Quickly disassembles into five, lightweight pieces for easy transport!

Wireless, quick disconnect in three easy steps:

Features

- Exclusive interlocking teeth for fully connectorless assembly/disassembly
- 250 lb. weight capacity
- Interchangeable color panels included with each unit—FREE!
- Auto grade cut pile carpet
- Tiller mounted basket

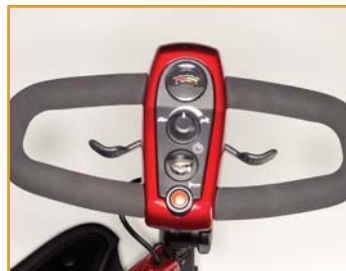

New control panel design and Delta tiller standard



Step 1: Pop off the seat!



Step 2: Pull off the battery pack!



Step 3: Push apart the frame!

* Battery range will vary due to rider weight, drive surface and drive slope.
** Electronics warranty excludes batteries.

Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.