

Med-Aire Plus Alternating Pressure Mattress Replacement System with Low Air Loss

14029 80" (L) x 36" (W) x 8" (H), Mattress and Pump System, 1/ea

HCPCS code: E0277

Available
1 Pack

BARIATRIC

450 LB

WEIGHT CAPACITY

FEATURES:

- This combination therapy system provides both alternating pressure and low air loss to optimize pressure redistribution and manage skin maceration.
- Designed to prevent, treat and manage pressure ulcers in the home, long term care, or acute care setting.
- **Cell-on-cell design** of 8" deep air cells prevents "bottoming out" and also provides up to 24 hours of power outage protection.
- **Digital pump** allows for customization to individual patients' comfort and safety needs, and will auto return to alternating pressure mode if system is left on static mode for >2 hours.
- **4-way stretch cover** is low shear, fluid-resistant, vapor permeable, quilted and zippered.
- Each one of the 20, **8" deep air bladders** are easily removed and replaced
- Pillow function maintains the 3 air cells at the head of the bed in static mode for patient comfort.
- 3 diagnostic alarms tell you whether a problem is with the mattress or the pump for more efficient servicing.
- 450 lb. weight capacity.
- CPR valve allows for rapid deflation.



14029



14427 (Electric Pump)

PARTS:

- 14029XPM 80" (L) x 36" (W) x 8" (H),
Mattress Only for 14029 with 3 Ports, 1/ea
- 14029XP 14029-CELL Pump Only for 14029, 1/ea (Shown)
- 14029-CELL Replacement Bladder for 14029, 1/ea
- 14029-CELL 1 Replacement Bladder with Laser Holes for 14029, 1/ea

Mattress:

- Dimensions: 80" (L) x 36" (W) x 8" (H)
- Base Material: zippered vinyl
- **Maximum Support Capacity: 450 lbs.**
- Limited 18 Month Warranty

SPECIFICATIONS:


Pump:

- Power: 120 VAC, 60 Hz
- Pressure Range: 20 mmHg to 60 mmHg (inner mattress pressure)
- Auto firm (60 mmHg, Auto return in 30 min.)
- Seat function (Pressure preset + 10 mmHg)
- Static selectable with auto-return to alternating mode.
- Adjustable cycle time 10, 15, 20, 25 minutes.
- Function lock (Auto lock in 30 seconds)
- Low pressure alarm
- System (alternating) alarm
- Power failure alarm
- Alarm mute
- Memory function
- Air Flow: 8-LPM
- Power Cord Length: 15 Ft.
- Pump Dimensions: 4.5" (L) x 14.2" (W) x 9.8" (H)
- UL: 60601-1, 1st Edition-2003
- Limited 18 Month Warranty

**NEW
DIGITAL
PUMP
DESIGN!**



14029XP

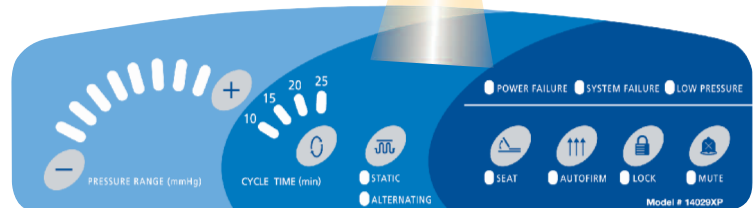

Important Notice: Item numbers that display an asterisk () reflect an item that will be discontinued or package change.