

# Med-Aire Plus Alternating Pressure Mattress Replacement System with Low Air Loss with Defined Perimeter

Available  
**1 Pack**

14029DP

80" (L) x 36" (W) x 10" (H), Pump and Mattress System, 1/ea

HCPSC code: **E0277**

BARIATRIC  
**450 LB**  
WEIGHT CAPACITY

## FEATURES:

- This bariatric system provides alternating pressure to optimize pressure redistribution for patients with up to Stage IV pressure ulcers in the home, long term care, or acute care setting.
- **10" deep static perimeter tubing surrounds the length of the mattress preventing the patient from becoming trapped between the mattress and side rails.**
- **Cell-on-cell design** of 8" deep air cells prevents "bottoming out" and also provides up to 24 hours of power outage protection.
- **Digital pump** allows for customization to individual patients' comfort and safety needs, and will auto-return to alternating pressure mode if system is left on static mode for >2 hours.
- **4-way stretch cover** is low shear, fluid-resistant, vapor permeable, quilted and zippered.
- Each one of the 20, **8" deep air bladders** are easily removed and replaced.
- Pillow function maintains the 3 air cells at the head of the bed in static mode for patient comfort.
- 3 diagnostic alarms tell you whether a problem is with the mattress or the pump for more efficient servicing.
- 450 lb. weight capacity.
- CPR valve allows for rapid deflation.



14029DP



14427 (Electric Pump)

## PARTS:

14029DPM

80" (L) x 36" (W) x 8" (H),  
Mattress Only, 1/ea

14029XP

Pump only, 1/ea

14029DP-CELL

Replacement Bladder for 14029DP, 1/ea

14029DP-CELL 1

Replacement Bladder with Laser Holes for 14029DP, 1/ea

## Mattress:

- Dimensions: 80" (L) x 36" (W) x 10" (H)
- 10" static elevated perimeter
- Base Material: zippered vinyl
- **Maximum Support Capacity: 450 lbs.**
- Limited 18 Month Warranty


## SPECIFICATIONS:

### Pump:

- Power: 120 VAC, 60 Hz
- Pressure Range: 20 mmHg to 60 mmHg (inner mattress pressure)
- Auto firm (60 mmHg, Auto return in 30 min.)
- Seat function (Pressure preset + 10 mmHg)
- Static selectable with auto-return to alternating mode.
- Adjustable cycle time 10, 15, 20, 25 minutes.
- Function lock (Auto lock in 30 seconds)
- Low pressure alarm
- System (alternating) alarm
- Power failure alarm
- Alarm mute
- Memory function
- Air Flow: 8-LPM
- Power Cord Length: 15 Ft.
- Pump Dimensions: 14.2" (W) x 9.8" (H) x 4.5" (L)
- UL: 60601-1, 1st Edition-2003
- Limited 18 Month Warranty



14029XP

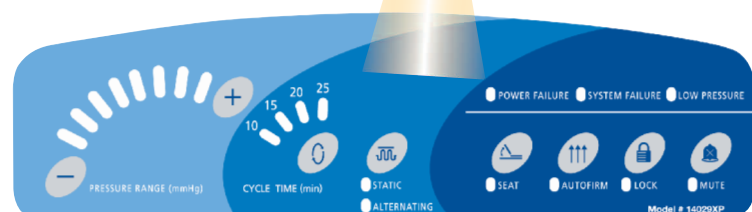

\*Important Notice: Item numbers that display an asterisk (\*) reflect an item that will be discontinued or package change.

To order call toll free: 877.224.0946

PRESSURE PREVENTION

drive™