


Gravity 8 Deluxe Long Term Care Pressure Redistribution Mattress


The Gravity 8 multi-layered foam system provides outstanding zoned pressure redistribution, with an emphasis on managing patients at risk for heel pressure ulcers. With a gently sloping heel section, and channel cuts in the middle layer of foam under the foot section, heel pressure off-loading and comfort are maximized.

NEW!

350 LB
WEIGHT CAPACITY

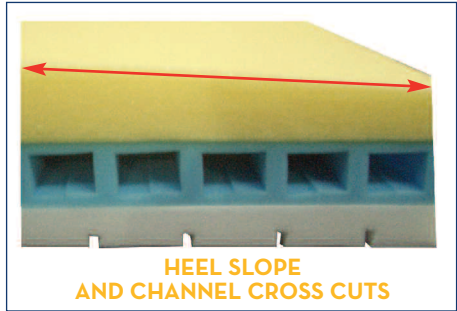

**PREVENT HEEL PRESSURE ULCERS
WITH THE GRAVITY 8'S HEEL SLOPE,
CHANNEL CUTS, AND OPEN-CELL
VISCO-ELASTIC MEMORY FOAM
FOOT SECTION!**

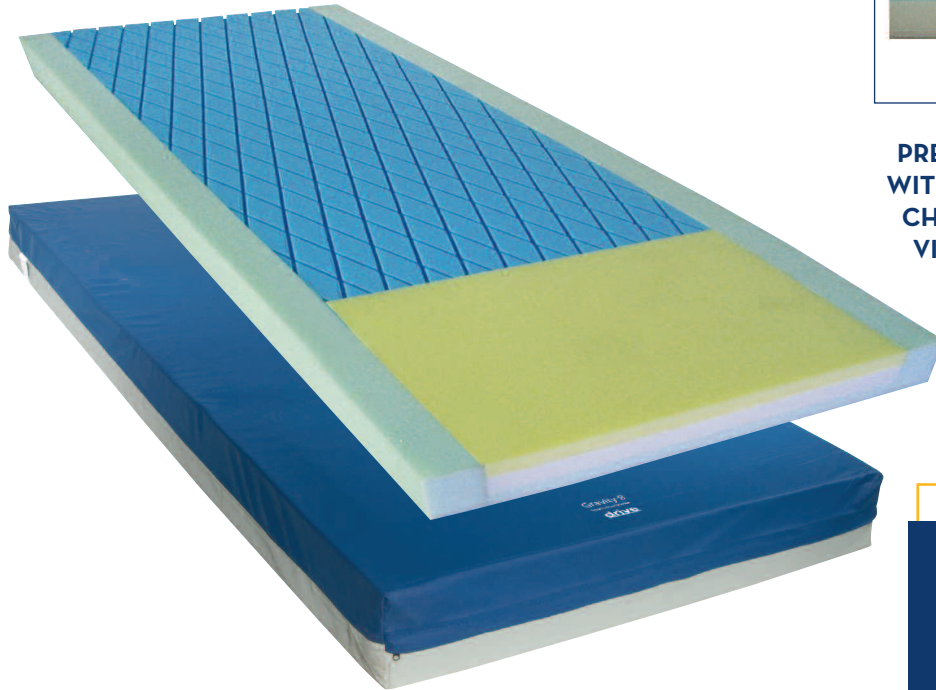

COMPARE TO:

Direct Supply® Panacea
Invacare® Therapy 2000
Mason® 6500
Span America® GeoMatt Max
Blue Chip® Elite

CUSTOM SIZES AVAILABLE

- 1P** 15870 80" (L) x 36" (W) x 6" (H), 1/ea
- 1P** 15888 80" (L) x 36" (W) x 6" (H) with Raised Side Rails, 1/ea
- 1P** 15884 84" (L) x 36" (W) x 6" (H), 1/ea
- 1P** 15885 84" (L) x 36" (W) x 6" (H) with Raised Side Rails, 1/ea
- 1P** 15886 76" (L) x 36" (W) x 6" (H), 1/ea
- 1P** 15887 76" (L) x 36" (W) x 6" (H) with Raised Side Rails, 1/ea

HCPCS code: E0184

FEATURES:

- Outstanding combination of multi-layered die-cut and memory foam available with or without raised side rails.
- ◆ Mattress includes a 6" to 5" heel slope and a "channel cut" center layer of foam under the foot section to target heel ulcer prevention.
- ◆ Top layer uses unique die-cutting and a true memory foam heel section to optimize zoned pressure redistribution.
- ◆ Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress significantly.
- ◆ Cover is fire retardant, fluid proof, tear resistant, low sheer, and easy to clean.
- ◆ **1 Year Warranty on the cover.**
- ◆ **Limited 10 Year Warranty.**

Important Notice: Item numbers that display an asterisk () reflect an item that will be discontinued or package change.

1P = 1-Pack Available