



RB14 *MEDITERRANEAN* Bariatric Side-Entry Bathtub

The RB14 Mediterranean is designed specifically for the larger person that needs extra room to bathe. The low threshold and wide side-entry door allows the user to enter and exit the bathtub easily and safely. The 28" extra-wide contoured seat allows the larger person to sit comfortably and safely in a slightly reclined position.

The sleek exterior design of the RB14 Mediterranean will give your bathroom an upscale look. Featuring a full range of standard and optional features, the RB14 Mediterranean is designed to enhance your bathing experience.

The Rane Bathing Systems product line provides a unique combination of style, comfort, ease-of-use, and low cost of ownership.

- ◆ **Full set of Standard Features**, including a Roman Deck Mount Valve Package made of solid brass with a polished chrome finish, Handheld Shower Wand, Removable Access Panel for use of transfer devices, and Side Panels for freestanding or corner installations.
- ◆ **A "Best in Class" Warranty** resulting in low cost of ownership.
- ◆ **Nationwide Distribution and Customer Service.**
- ◆ **Made in the USA.** Manufactured entirely in Sparta, TN from the highest quality components and materials.



MADE IN THE USA



"NOTHING IS MORE REFRESHING THAN RANE"

Standard Features


- ◆ **Roman Deck Mount Valve Package** made of solid brass with a polished chrome finish that provides an elegant look to your bathroom.
- ◆ **Handheld Shower Wand** for hair and upper body washing.
- ◆ **Removable Access Panel** for use of transfer devices.
- ◆ **Side Panels** for both the left or right side of the tub, which is needed for freestanding or corner installations.



The 28" extra-wide contoured seat is designed to allow persons of size the opportunity to bathe in comfort.

Optional Features

- ◆ **CleanRane™ Air Spa System** with variable speeds that enhance the therapeutic and relaxing effects of bathing for the user.
- ◆ **Supreme Combination Package** featuring CleanRane™ Air Spa System, Heated Whirlpool System, Aromatherapy, and Chromatherapy.
- ◆ **Brushed Nickel or Oil Rubbed Bronze Roman Deck Mount Package.**
- ◆ **3/4" Supply Lines** constructed of stainless steel mesh clad.
- ◆ **Head Pillow** with suction cups for comfort and head support.
- ◆ **Gel Seat Cushion** made with an ultra-soft viscoelastic polymer gel and custom designed to fit perfectly in our residential tubs.
- ◆ **Lap and Shoulder Harness** for user safety.
- ◆ **Secondary Door Lock** to prevent accidental door opening when in use.
- ◆ **Chrome Shower Bar** to install shower wand on wall, which will allow stand-up showering if desired.

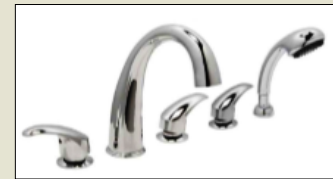

Roman Deck Mount Valve Package made of solid brass with a polished chrome finish. Available in Brushed Nickel or Oil Rubbed Bronze as an upgrade.

Why CleanRane™ Air Spa System?

The CleanRane™ Air Spa System enhances the therapeutic and relaxing effects of bathing and provides a pleasant experience for the resident. Hydromassage provides millions of preheated air bubbles in the water which increases stimulation of peripheral blood circulation, reduces pain and inflammation in the joints, gently massages away dirt, and induces deeper relaxation. The motor is quiet, efficient, and has a ceramic heating element that blows warm air to help maintain the water temperature and warm the surface of the tub.



The CleanRane™ Air Spa System features 12 jets strategically placed for maximum effect.

Technical Specifications

Width:	33.5"	Weight capacity (occupant):	600 lbs. / 272.7 kg
Length:	60"	Water capacity (maximum):	68 gallons / 257.4 L
Height:	36.5" (lowest setting)	Fill time (68 gallons):	5 minutes, 40 seconds
Door Width:	41"	Drain time (68 gallons):	2 minutes, 15 seconds
Adjustable leveling legs:	2" maximum	Water supply inlet:	3/4" diameter
Seat width (interior):	28"	Drain size:	1.5"
Seat height:	21"	Electrical supply:	115V AC; 60Hz; 15 amp
Threshold height:	13" (lowest setting)	Door configuration:	Outward swing, Right or Left Hinge
Transfer Device Access:	7" H x 41" W		

Certifications: Plumbing; IAPMO Z124-2011, CSA B45.5-11 / Electrical; UL 1795, CSA STD C22.2 NO. 218.2

Warranty: **5 yrs.** for external finishes and composite shell, **3 yrs.** for all moving and stationary parts, and **lifetime** on the door seal.



Rane Bathing Systems | 155 Norcom Way, Sparta, TN 38583 | Phone: 888.880.7373 | Fax: 931.836.6374

www.Ranetubs.com