### **New ADA Laws**

- All of our lifts are compliant with the New ADA Laws
- All swimming pools in the USA that fall under Title II or Title III need to comply with the new laws. Title II (Public Industry) municipal pools, school pools,

government owned pools, etc.


Title III (Private Industry) place of recreation, place of lodging, which addresses public accommodations

- Existing pools shall become compliant with accessibility by March 15, 2012
- There are tax incentives available See details at www.GlobalLiftCorp.com

PLEASE CHOOSE WHICH APPLICATION YOUR POOL GUTTER BEST REPRESENTS




### **DECK CONFIGURATOR COMMERCIAL SERIES LIFTS**

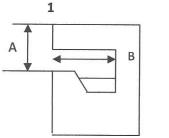


**Fully Recessed Gutter** 




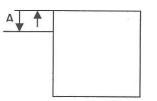

Gutter less Deck

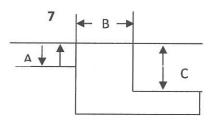

Above Ground Pool

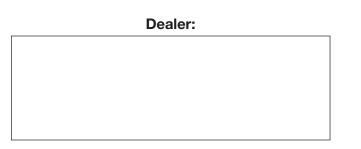

NOTE: If your pool doesn't fall within these measurements,

Please contact your dealer for assistance

**LEGEND MIN/MAX** A= (Draft) pool deck to water line 0 – 7"

B= Width of Gutter 1 – 12" C= Height of Curb 1 – 2"

## **Commercial Series Lifts**



### **Commercial Series Lifts**

#### **C-375**

#### C-450

### **Standard Features**

- ADA Compliant Pool Lift
- Stainless Steel Construction
- Powder Coat Finish
- Lift Capacity -C-375 375 LBS. / C-450 450 LBS.
- Battery Powered 24 volt system
- User Friendly controls
- Lifetime Structural Warranty
- Arm Rests
- Foot Rest
- Safety Seat Belt
- Right Seat Position
- UPS Packaging for easy shipping
- Can be used for spas and pools
- Fits Most Standard Pool Applications
- Key Lock Out System (Patent Pending)
- Emergency Return System (Patent Pending)
- UV Resistant
- Includes Battery Charging Station
- Retro fit, or Tri-Point drop in, or Tri-Point Anchoring System optional

#### **Global Lift Corp. Products are** Laboratory tested to meet ADA **Standards**





### www.GlobalLiftCorp.com

# **Commercial Series Lifts**

- eliminates pinch points.
- portable unit with the addition of a portable kit at any time.
- turn of a key when the unit is not in or intended to be in use.
- out of the water.

**Options** 



The new Anti-Pinch Linkage System (Patent Pending) holds the seat in a level position with an external linkage system that virtually

Portable Capability System (Patent Pending) – The commercial series lifts can be purchased as a stationary lift and can be made into a fully

The key lock out system is designed to lock the system out with the

The emergency return system is designed to connect the power source directly to the actuator in the event the handset and or controls malfunction. It will safely return the stranded user to the up position and



**Extra Battery** 



**Commercial Series Portable Kit** with Quick Attach Bolts