

Hoyer® HPL700

The HPL700 is a uniquely designed electronically operated patient lift for the safe lifting of heavier patients. It is both sturdy and robust with a safe working load rating of 700 lbs. The HPL700 has been designed to be configured with or without a weigh scale. The weigh scale is seamlessly integrated with a digital display on both sides and user-friendly controls.




6-Point Cradle

The large, open cradle results in less discomfort to the patient caused by pinching and squeezing.



Powered Base

The powered base allows the legs to be opened and closed to enable access around armchairs, wheelchairs and similar obstacles without bending or twisting.

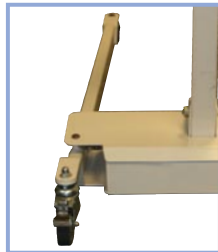

Power Source

The HPL700 comes standard with two batteries and a charging station. Visual and audible low-battery indicators are featured.



Optional Scale

The optional scale allows the caregiver to adjust the position of the display for easier viewing. Simple, intuitive controls provide the resident's weight in pounds or kilograms.



Specification	Imperial	Metric
Safe Working Load	700 lbs	227kg
Maximum Overall Length	56"	1423mm
Minimum Overall Length	55"	1397mm
Maximum Overall Height	81"	2058mm
Minimum Overall Height	57.75"	1467mm
Height at Maximum Reach	46.5"	1181mm
Reach at Maximum Height	28.5"	724mm
Reach at Minimum Height	27"	686mm
Maximum Reach*	37.75"	959mm
Turning Radius	62"	1575mm
Legs Open - External Width	46.25"	1175mm
Legs Open - Internal Width	42.25"	1073mm
Legs Closed - External Width	27.5"	699mm
Legs Closed - Internal Width	23.5"	597mm
Overall Height of Legs	4.5"	114mm
Ground Clearance	0.625"	16mm
Front Twin Casters	3"	76mm
Rear Braked Casters	4"	102mm

*Reach = Center of spreader bar to the front of the mast

Weights	Imperial	Metric
Mast, Base, Boom Assembly, Cradle	143.4 lbs	65kg
Power Pack	6.6 lbs	3kg
Total	150 lbs	68kg
Base Assembly	83 lbs	38kg
Mast and Boom (not including battery)	40.4 lbs	18kg
6-Point Spreader Bar	20 lbs	9kg

All measurements are within a +5/-5% tolerance.

6-Point Spreader Bar (measurement to top of location pins)

Cradle Maximum Height	69"	1753mm
Cradle Minimum Height	23.25"	591mm

Electric Shock Protection


Charger – Class II
Lift – Internal power source

Degree of Shock Protection

Charger – Type B
Lift – Type B
Intended operating environment: >+5°C <+40°C
Outside this environment, functionality and safety may be compromised

Electrical Specifications

Batteries – 2x12 volt Rechargeable sealed lead acid type
Battery Capacity – 3.2 Ampere Hours
Charger Rated Input – 100/240VAC 50/60Hz
Charger Rated Output – 27.4/29 VDC@0.8A

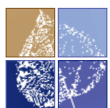

Joerns
Healthcare

www.joerns.com • email: info@joerns.com

Post Acute, Acute, HomeCare
5001 Joerns Drive
Stevens Point WI 54481
(P) 800.826.0270
(F) 800.457.8827

VA/Government
7027 Hayvenhurst Avenue
Van Nuys, CA 91406
(P) 800.966.6662
(F) 800.232.9796

Canadian Office
1000 Clarke Road, Suite 6
London, ON Canada N5V 3A9
(P) 866.546.1151
(F) 519.451.8662