

Joerns Healthcare is committed to providing a complete line of top quality lifting and repositioning equipment to reduce workplace injury risks and overall facility insurance costs. Our Hoyerpro+ lifts have been designed to work together to *Create a Safer Environment* (C.A.S.E.) for your staff and patients, while addressing your budget requirements.

The Hoyerpro+ 1000 Ceiling Lift creates a superior solution for your bariatric lifting and repositioning needs. It features an industry leading lifting capacity of up to 1000 lbs. to allow for the safe transfer of your heaviest patients. A high output lifting motor allows the Hoyerpro+ 1000 to perform all lifts with one system.


Heavy-duty motors and all-metal drive trains ensure long term reliability and safety. The high capacity batteries allow for more lifts per charge than regular ceiling lifts. An extra wide, extra long lifting strap provides for maximum safety and convenience; floor lifts can easily be accomplished in rooms with high ceilings. The standard oversize spreader bar helps provide maximum patient comfort.

Joerns® Hoyerpro+ 1000 Ceiling Lift Features

- ① Industry exclusive trolley design accommodates curtain track solution.
- ② Up to 1,000lb lifting capacity. Choose manual or power traverse.
- ③ Compact design with high output motor design.
- ④ Integrated safety features include slack strap cutout switch, angle detector, anti-freewheel, limit switches, overload protection, emergency mechanical lowering and integrated cycle counter.
- ⑤ Extra wide and extra long lifting strap for safety and convenience.
- ⑥ Super six point padded spreader bar for superior safety and comfort.
- ⑦ Emergency on/off.
- ⑧ 2 or 4 function hand controller electric operation for optimal reliability.
- ⑨ Smart sleep mode function to conserve battery power.
- ⑩ Digital control pad with visual and audio low battery indicators.
- ⑪ All metal lifting and traverse drive train for safety, reliability and durability.

Options

Digital Weigh Scale

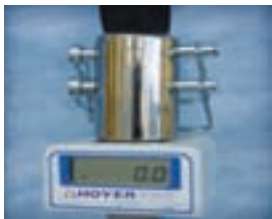

The optional digital weigh scale allows for accurate patient weighing in lbs. and kgs.

Curtain Track Pass Through



A patent pending trolley and track design allows your existing privacy curtain tracks to work seamlessly with ceiling lift tracks.

You must use Hoyer® slings with Hoyer lifts.
Specifications may vary due to manufacturing tolerances

Specifications

| | |
|--------------------------------|-------------------------------------|
| Power Source: | 2x12V SLA Batteries, 7.2 AH |
| Safe Working Load | |
| Model 1000/1000M: | 1000 lbs (454 kgs) |
| Height (excluding trolley): | 10.5" (26.2 cm) |
| Width: | 14.25" (35.6 cm) |
| Length: | 18" (45 cm) |
| Spreader Bar: | 6 point with auto leveling sidebars |
| Lifting Strap Width: | 2" (5 cm) |
| Lifting Strap Length: | 8' (243 cm) |
| Charger: | 120VAC/24VDC CSA/UL Approved |
| Lifting Motor: | 24VDC |
| Lifting Motor Drive System: | Direct Drive Gear |
| Traverse Motor (T7003): | 24 VDC |
| Traverse Motor Drive System: | Chain/Sprocket |
| Overall Weight: | 47 lbs (21.4 kgs) |
| Charging: | End of Track |
| Emergency Stop: | Standard |
| Visual Battery Indicator: | Standard |
| Audio Low battery Indicator: | Standard |
| Emergency Lowering: | Standard |
| Slack Strap Detector: | Standard |
| Angle Detector Cutout Switch: | Standard |
| Anti Free Wheel Safety System: | Standard |
| Up Limit Switch: | Standard |
| Down Limit Switch: | Standard |
| Remote hand Controller: | Standard |
| Warranty: | 1 Year |
| Digital Weigh Scale: | Optional |
| Approved: | CSA-Z10535-03 |



You must use Hoyer® slings with Hoyer lifts.
Specifications may vary due to manufacturing tolerances

© 2008 Joerns Healthcare • Part# 6190081 RevA • DCO# 08-0733

Joerns Healthcare Inc.
5001 Joerns Drive, Stevens Point WI 54481
(P) 800.826.0270 • (F) 800.457.8827