

DEALER			COLOURS INFORMATION For Office Use Only		
Dealer:			Order Date:		
Account: P.O.#:			S/N:		
Ordered By:			RTS		
Address:			RTS: Phone:		
City: State: Zip:			Email:		
Phone: Fax:			SHIP TO		
Email:			Name: Ship Via:		
Customer Name:			Address:		
Customer Weight:			City: State: Zip:		
			Phone: Fax:		

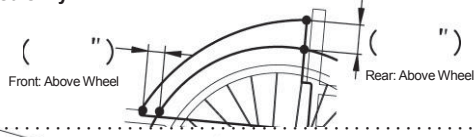
CHAIR INFORMATION	PRICE
CHUMP E1235	Base Wheel Chair Price E1235 \$ 2,400
Adjustable FSH-RSH	Total Section A \$
Max. Capacity 250Lbs	Total Section B \$
NOTE	Total Section C \$
Main Frame : ZK55 & ZK60 High Grade A7000 Series Aluminum	Total Section D \$
Main Frame : Lifetime Warranty	Total Section E \$
*See warranty policy	Total Section F (Custom Parts Color) \$
	ORDER GRAND TOTAL \$

A FRAME COLORS & PARTS COLORS

FRAME COLORS		
Powder-Co ☐ Gloss Black ☐ Twilight Flake ☐ Yellow ☐ Magenta ☐ Fire Engine Red ☐ White ☐ Candy Blue ☐ Candy Purple ☐ Candy Green ☐ Candy Red ☐ Candy Raspberry ☐ Cobalt Blue ☐ Matte Gray ☐ Matte Purple ☐ Matte Blue ☐ Matte Green ☐ Matte Black ☐ Matte Red ☐ Neon Red ☐ Neon Green ☐ Orange Sparkle ☐ Mellon Yellow ☐ Gun Smoke Silver ☐ Pigeon Blue ☐ Granny Sparkle Green ☐ Sparkle Sky Blue ☐ Vermilion ☐ Reed Green ☐ Blue Lilac		\$ STD
Special ☐ Polished ☐ 24 karat Gold Plated		\$ 800 \$ CFP
PARTS COLORS (Custom Parts Color Section, Please see Section F)		
STD Color	Black Color (*Suspension Springs-Silver Color)	\$ STD
Custom Color Main Bolt on Parts One Color Package Including: Backrest, Foot Plate, Footrest Rings, Backrest Brackets, Saddle Stops, Forks & Knurl Nuts, Caps, Front Suspension Springs, Axle Tube, Axle Clamps, Axle Plates, Axle Bases, Camber Discs, Caster Housings ☆marked parts	☐ Powder-Coated () ☐ Polished ☐ 24 karat Gold Plated	\$ 500 \$ 750 \$ CFP
Total Section A		\$

Chump: Page 2 of 9

B FRAME MEASUREMENTS 2/2

5. BACKREST HEIGHT * Rigidizer bar 3.25" deep std.		Adj. Range		Adj. Range		Adj. Range				\$ STD																																																																																																																
		☐ 8"	☐ 9"	☐ 10"	☐ 11"	☐ 12"	☐ 13"	☐ 14"	☐ 15"	☐ 16"	☐ 17"	☐ 18"																																																																																																														
Type		☐ Locking mechan. with adjust. height & angle \$ STD ☐ Double Locking mechan. with adjust. height & angle \$ 25 ☐ Folding with adjust. height & angle \$ NCO																																																																																																																								
CUSTOM		☐ Welded Backrest Angle Set at (") to the ground *Non Folding Backrest Only \$ 150 ☐ Welded Side Guards *Non Folding Backrest Only \$ 450  ☐ Custom Backrest \$ 150 (ie: double rigidizer bars, custom depth, custom height, etc)																																																																																																																								
6. BACKREST ANGLE		Set at (°) to the ground *STD 90°									\$ STD																																																																																																															
7. FRONT SEAT HEIGHT		☐ 14.5" ☐ 15.0" ☐ 15.5" ☐ 16.0" ☐ 16.5" ☐ 17.0" ☐ 17.5" ☐ 18.0" ☐ 18.5" \$ STD The front seat heights are impacted by the size of your front caster.																																																																																																																								
Front Seat Height Ranges *=Not Available with Quick Release Caster Axle.		<table border="1"> <thead> <tr> <th colspan="2" rowspan="2"></th> <th colspan="9">Front Seat Height</th> </tr> <tr> <th>14.5"</th> <th>15.0"</th> <th>15.5"</th> <th>16.0"</th> <th>16.5"</th> <th>17.0"</th> <th>17.5"</th> <th>18.0"</th> <th>18.5"</th> </tr> </thead> <tbody> <tr> <td rowspan="8">Your Casters & Forks</td> <td>3"Caster & STD Forks</td> <td>☐ *</td> <td>☐ *</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> </tr> <tr> <td>3"Caster & Springs</td> <td></td> <td></td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> </tr> <tr> <td>3"Caster & Frog Legs</td> <td></td> <td></td> <td></td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> </tr> <tr> <td>4"Caster & STD Forks</td> <td></td> <td></td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> </tr> <tr> <td>4"Caster & Springs</td> <td></td> <td></td> <td></td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> </tr> <tr> <td>4"Caster & Frog Legs</td> <td></td> <td></td> <td></td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> </tr> <tr> <td>5"Caster & STD Forks</td> <td></td> <td></td> <td></td> <td></td> <td>☐ *</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐ *</td> </tr> <tr> <td>5"Caster & Springs</td> <td></td> <td></td> <td></td> <td></td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> </tr> <tr> <td>5"Caster & Frog Legs</td> <td></td> <td></td> <td></td> <td></td> <td></td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> <td>☐ *</td> </tr> </tbody> </table>												Front Seat Height									14.5"	15.0"	15.5"	16.0"	16.5"	17.0"	17.5"	18.0"	18.5"	Your Casters & Forks	3"Caster & STD Forks	☐ *	☐ *	☐	☐	☐	☐	☐	☐	☐	3"Caster & Springs			☐	☐	☐	☐	☐	☐	☐	3"Caster & Frog Legs				☐ *	☐ *	☐ *	☐ *	☐ *	☐ *	4"Caster & STD Forks			☐	☐	☐	☐	☐	☐	☐	4"Caster & Springs				☐	☐	☐	☐	☐	☐	4"Caster & Frog Legs				☐ *	☐ *	☐ *	☐ *	☐ *	☐ *	5"Caster & STD Forks					☐ *	☐	☐	☐	☐ *	5"Caster & Springs					☐ *	☐ *	☐ *	☐ *	☐ *	5"Caster & Frog Legs						☐ *	☐ *	☐ *	☐ *
		Front Seat Height																																																																																																																								
		14.5"	15.0"	15.5"	16.0"	16.5"	17.0"	17.5"	18.0"	18.5"																																																																																																																
Your Casters & Forks	3"Caster & STD Forks	☐ *	☐ *	☐	☐	☐	☐	☐	☐	☐																																																																																																																
	3"Caster & Springs			☐	☐	☐	☐	☐	☐	☐																																																																																																																
	3"Caster & Frog Legs				☐ *	☐ *	☐ *	☐ *	☐ *	☐ *																																																																																																																
	4"Caster & STD Forks			☐	☐	☐	☐	☐	☐	☐																																																																																																																
	4"Caster & Springs				☐	☐	☐	☐	☐	☐																																																																																																																
	4"Caster & Frog Legs				☐ *	☐ *	☐ *	☐ *	☐ *	☐ *																																																																																																																
	5"Caster & STD Forks					☐ *	☐	☐	☐	☐ *																																																																																																																
	5"Caster & Springs					☐ *	☐ *	☐ *	☐ *	☐ *																																																																																																																
5"Caster & Frog Legs						☐ *	☐ *	☐ *	☐ *																																																																																																																	
8. REAR SEAT HEIGHT		☐ 12.5" ☐ 13.0" ☐ 13.5" ☐ 14.0" ☐ 14.5" ☐ 15.0" ☐ 15.5" ☐ 16.0" ☐ 16.5" ☐ 17.0" ☐ 17.5" \$ STD																																																																																																																								
Rear Seat Height Ranges		<table border="1"> <thead> <tr> <th>Rear Tire Size</th> <th>Rear Seat Height</th> </tr> </thead> <tbody> <tr> <td>20"x1"(451)</td> <td>12.5"- 16.0"</td> </tr> <tr> <td>20"x1 3/8"(451)</td> <td>13.0"- 16.5"</td> </tr> <tr> <td>22"x1"(501)</td> <td>13.5"- 17.0"</td> </tr> <tr> <td>22"x1 3/8"(501)</td> <td>14.0"- 17.5"</td> </tr> </tbody> </table>										Rear Tire Size	Rear Seat Height	20"x1"(451)	12.5"- 16.0"	20"x1 3/8"(451)	13.0"- 16.5"	22"x1"(501)	13.5"- 17.0"	22"x1 3/8"(501)	14.0"- 17.5"																																																																																																					
Rear Tire Size	Rear Seat Height																																																																																																																									
20"x1"(451)	12.5"- 16.0"																																																																																																																									
20"x1 3/8"(451)	13.0"- 16.5"																																																																																																																									
22"x1"(501)	13.5"- 17.0"																																																																																																																									
22"x1 3/8"(501)	14.0"- 17.5"																																																																																																																									
9. CENTER OF GRAVITY		☐ 0.0" ☐ 0.5" ☐ 1.0" ☐ 1.5" ☐ 2.0" ☐ 2.5" ☐ 3.0" ☐ 3.5" ☐ 4.0" ☐ 4.5" \$ STD																																																																																																																								
Center of Gravity Ranges ☆=Not Available 4" and 5" Caster. ★=Not Available 22" Wheel.		<table border="1"> <thead> <tr> <th colspan="2" rowspan="2"></th> <th colspan="9">Your Center of Gravity</th> </tr> <tr> <th>0"</th> <th>0.5"</th> <th>1.0"</th> <th>1.5"</th> <th>2.0"</th> <th>2.5"</th> <th>3.0"</th> <th>3.5"</th> <th>4.0"</th> <th>4.5"</th> </tr> </thead> <tbody> <tr> <td rowspan="5">Your Seat Depth</td> <td>10"</td> <td>☐ ☆</td> <td>☐ ☆★</td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> </tr> <tr> <td>11"</td> <td>☐ ☆</td> <td>☐ ☆</td> <td>☐ ☆</td> <td>☐ ☆★</td> <td></td> <td></td> <td></td> <td></td> <td></td> <td></td> </tr> <tr> <td>12"</td> <td>☐</td> <td>☐</td> <td>☐ ☆</td> <td>☐ ☆</td> <td>☐ ☆</td> <td>☐ ☆★</td> <td></td> <td></td> <td></td> <td></td> </tr> <tr> <td>13"</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐ ☆</td> <td>☐ ☆</td> <td>☐ ☆</td> <td>☐ ☆★</td> <td></td> <td></td> </tr> <tr> <td>14"</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐ ☆</td> <td>☐ ☆</td> <td>☐ ☆★</td> </tr> <tr> <td></td> <td>15"</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐</td> <td>☐ ☆</td> </tr> </tbody> </table>												Your Center of Gravity									0"	0.5"	1.0"	1.5"	2.0"	2.5"	3.0"	3.5"	4.0"	4.5"	Your Seat Depth	10"	☐ ☆	☐ ☆★									11"	☐ ☆	☐ ☆	☐ ☆	☐ ☆★							12"	☐	☐	☐ ☆	☐ ☆	☐ ☆	☐ ☆★					13"	☐	☐	☐	☐	☐ ☆	☐ ☆	☐ ☆	☐ ☆★			14"	☐	☐	☐	☐	☐	☐	☐	☐ ☆	☐ ☆	☐ ☆★		15"	☐	☐	☐	☐	☐	☐	☐	☐	☐	☐ ☆																						
		Your Center of Gravity																																																																																																																								
		0"	0.5"	1.0"	1.5"	2.0"	2.5"	3.0"	3.5"	4.0"	4.5"																																																																																																															
Your Seat Depth	10"	☐ ☆	☐ ☆★																																																																																																																							
	11"	☐ ☆	☐ ☆	☐ ☆	☐ ☆★																																																																																																																					
	12"	☐	☐	☐ ☆	☐ ☆	☐ ☆	☐ ☆★																																																																																																																			
	13"	☐	☐	☐	☐	☐ ☆	☐ ☆	☐ ☆	☐ ☆★																																																																																																																	
	14"	☐	☐	☐	☐	☐	☐	☐	☐ ☆	☐ ☆	☐ ☆★																																																																																																															
	15"	☐	☐	☐	☐	☐	☐	☐	☐	☐	☐ ☆																																																																																																															
10. CAMBER ANGLE		☐ 2° ☐ 4° ☐ 6° \$ STD																																																																																																																								
11. FRONT OF SEAT SLING TO FOOTREST		☐ 12.0" ☐ 12.5" ☐ 13.0" ☐ 13.5" ☐ 14.0" ☐ 14.5" ☐ 15.0" ☐ 15.5" ☐ 16.0" ☐ 16.5" ☐ 17.0" \$ STD *Measurement must be 1.5" shorter than your Front Seat Height																																																																																																																								
HIGH MOUNT ☐ 6.0" ☐ 6.5" ☐ 7.0" ☐ 7.5" ☐ 8.0" ☐ 8.5" ☐ 9.0" ☐ 9.5" ☐ 10.0" ☐ 10.5" ☐ 11.0" ☐ 11.5" \$ 65 With Risers or Clamps																																																																																																																										
12. INSIDE WHEEL CLEARANCE		☐ 3/4" * ☐ 1" ☐ 1 1/4" ☐ 1 1/2" ☐ 1 3/4" ☐ 2" \$ STD *3/4" Inside wheel clearance not available with 2 degree camber. (with the exception of Y-360 wheels) *3/4" Inside wheel clearance not available with Armrest *3/4" Inside wheel clearance not available with D's Locks																																																																																																																								
NOTE Chump is a custom made wheelchair. If you would like to customize your Chump beyond what is available on this order form please submit your request here. Colours will determine as to whether this request can be accomplished while maintaining the safety and integrity of the wheelchair. If it can be accomplished the quote will reflect any and all custom charges, thank you.																																																																																																																										
Sub Total Section B 1/2										\$																																																																																																																
Sub Total Section B 2/2										\$																																																																																																																
Total Section B										\$																																																																																																																

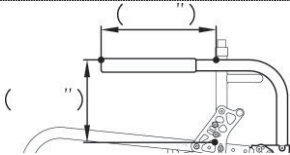
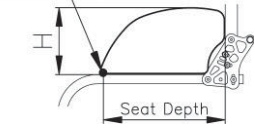
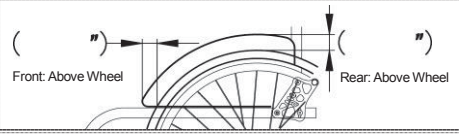
C WHEELS, TIRES & CASTERS, FORKS 1/2

WHEELS														
Standard Radial Spoke Wheel ☐ 20" (451) ☐ 22" (501) 		\$ STD												
Mag Wheel ☐ 20" (451) ☐ 22" (501) 		\$ STD												
w/Chrome Handrims (Long Tab Only)														
Twister Super Strength wheel ☐ 20" (451) ☐ 22" (501) 		\$ 300												
Spinergy Spox 18 Spokes	☐ 22" (501) Spoke Color ☐ Black ☐ Red ☐ Yellow ☐ Blue ☐ White 	\$ 695												
	With Silver Hub With Black Hub With Black Hub With Black Hub With Black Hub													
	☐ Custom Color Spoke Pattern <table><tr><th colspan="2">Left Wheel</th><th>Right Wheel</th></tr><tr><td>Spokes</td><td>Inside</td><td></td></tr><tr><td>Color</td><td>Outside</td><td></td></tr><tr><td>Hub Color</td><td colspan="2">☐ Black/☐ Silver ☐ Black/☐ Silver</td></tr></table>	Left Wheel		Right Wheel	Spokes	Inside		Color	Outside		Hub Color	☐ Black/ ☐ Silver ☐ Black/ ☐ Silver		\$ 795
	Left Wheel		Right Wheel											
	Spokes	Inside												
Color	Outside													
Hub Color	☐ Black/ ☐ Silver ☐ Black/ ☐ Silver													
Spinergy LX 12 Spokes	☐ 22" (501) Spoke Color ☐ Black ☐ Red ☐ Yellow ☐ Blue 	\$ 895												
	With Silver Hub With Black Hub With Black Hub With Black Hub With Black Hub													
	☐ Custom Color Spoke Pattern <table><tr><th colspan="2">Left Wheel</th><th>Right Wheel</th></tr><tr><td>Spokes</td><td>Inside</td><td></td></tr><tr><td>Color</td><td>Outside</td><td></td></tr><tr><td>Hub Color</td><td colspan="2">☐ Black/☐ Silver ☐ Black/☐ Silver</td></tr></table>	Left Wheel		Right Wheel	Spokes	Inside		Color	Outside		Hub Color	☐ Black/ ☐ Silver ☐ Black/ ☐ Silver		\$ 995
	Left Wheel		Right Wheel											
	Spokes	Inside												
Color	Outside													
Hub Color	☐ Black/ ☐ Silver ☐ Black/ ☐ Silver													
Pro Wheels ☐ 20" (451) ☐ 22" (501) 		\$ 825												
w/ Aluminum Color Anodized Handrim Spoke & Handrim Color ☐ Anodized Black ☐ Anodized Red ☐ Anodized Blue														
*Not available with D's Lock														
Wijit Drive System *Call for availability		\$ 5,488												
HANDRIMS														
Type	☐ Standard For 20", 22" Wheel Anodized Black *Not available with Mag wheel, Pro wheels 	\$ STD												
	☐ Vinyl Coated For 20", 22" Wheel *Not available with Mag wheel, Pro wheels 	\$ 142												
	☐ Natural Fit For 22" Wheel Only *not available with Mag wheel.													
	☐ Standard Grip (Not available with short tabs.) 	\$ 295												
	☐ Super Grip (Not available with short tabs.) 	\$ 340												
	☐ Without Thumb Grip 	\$ 270												
	☐ Natural Fit (LT) For 22" Wheel Only *not available with Mag wheel.													
	☐ Standard Grip (Not available with short tabs.) 	\$ 325												
	☐ Super Grip (Not available with short tabs.) 	\$ 370												
	☐ Without Thumb Grip 	\$ 300												
☐ The Surge Handrims *not available with Mag wheel. 	\$ 315													
☐ The Surge Handrims (LT) *not available with Mag wheel. 	\$ 340													
TAB LENGTH	☐ Long Tabs(STD inside clearance 3/4") 	\$ STD												
	☐ Short Tabs(Inside clearance 1/2") *Not available with Mag wheel 	\$ STD												
	☐ Short Tab (Set closest to wheel cut off excess) 	\$ 15												
	*Not available with Mag Wheel, Pro Wheel.													
NOTE														

C WHEELS, TIRES & CASTERS, FORKS 2/2

TIRES					
☐ Low Pressure Tire ☐ 20"x1-3/8" or ☐ 22"x1-3/8" Gray Only 					E2211/E2212 \$ STD
☐ High Pressure Tire 22"x1" Gray Tire with black wall Only 					E2211/E2212 \$ NCO
☐ Solid Tire ☐ 20"x1", ☐ 22"x1" ☐ Dark Gray ☐ Black(\$30up charge) 					E2218 \$ 135
☐ E2213 Flat-Free Inserts For 20"x1-3/8" or 22"x1-3/8" Tires Only 					E2211/E2213 \$ 117
CASTERS					
SIZE	TIRE WIDTH	HUB COLOR	TIRE COLOR	DESCRIPTION	
☐ 3"	1"	Black Hub	Black Tire Only	Polyeurathane	\$ STD
☐ 3"	1"	Plastic Black Hub	Clear Tire Only	Lighted, Polyeurathane	\$ 65
☐ 3"	1"	Red Hub	Red Tire Only	Polyeurathane	\$ STD
☐ 3"	1"	Black Hub	Clear Tire Only	Polyeurathane	\$ STD
☐ 3"	1"	Aluminum Silver Hub	Black Tire Only	Made By Flog Legs	\$ 165
☐ 3"	1"	Aluminum Black Hub	Black Tire Only	Made By Flog Legs	\$ 165
☐ 3"	1.5"	Aluminum Silver Hub	Black Tire Only	Soft Roll Made By Flog Legs	\$ 165
☐ 3"	1.5"	Aluminum Black Hub	Black Tire Only	Soft Roll Made By Flog Legs	\$ 165
☐ 4"	1"	Plastic Black Hub	Black Tire Only		\$ STD
☐ 4"	1"	Plastic Black Hub	Clear Only	Lighted, Polyeurathane	\$ 65
☐ 4"	1"	Plastic Black Hub	Black Tire Only	Made By Flog Legs	\$ 125
☐ 4"	1"	Aluminum Silver Hub	Black Tire Only	Made By Flog Legs	\$ 165
☐ 4"	1"	Aluminum Black Hub	Black Tire Only	Made By Flog Legs	\$ 165
☐ 4"	1.5"	Composite Black Hub	Tire ☐ Black ☐ Blue ☐ Red	Soft Roll Made By Flog Legs	\$ 125
☐ 4"	1.5"	Aluminum Silver Hub	Tire ☐ Black ☐ Blue ☐ Red	Soft Roll Made By Flog Legs	\$ 165
☐ 4"	1.5"	Aluminum Black Hub	Black Tire Only	Soft Roll Made By Flog Legs	\$ 165
☐ 5"	1"	Composite Black Hub	Black Tire Only		\$ STD
☐ 5"	1"	Plastic Black Hub	Clear Tire Only	Lighted, Polyeurathane	\$ 65
☐ 5"	1"	Composite Black Hub	Black Tire Only	Made By Flog Legs	\$ 125
☐ 5"	1"	Aluminum Silver Hub	Black Tire Only	Made By Flog Legs	\$ 165
☐ 5"	1"	Aluminum Black Hub	Black Tire Only	Made By Flog Legs	\$ 165
☐ 5"	1.5"	Composite Black Hub	Tire ☐ Black ☐ Blue ☐ Red	Soft Roll Made By Flog Legs	\$ 135
☐ 5"	1.5"	Aluminum Silver Hub	Tire ☐ Black ☐ Blue ☐ Red	Soft Roll Made By Flog Legs	\$ 165
☐ 5"	1.5"	Aluminum Black Hub	Black Tire Only	Soft Roll Made By Flog Legs	\$ 165
FORKS					
☐ STD Aluminum Forks Black Only 					\$ STD
Extra Options for Aluminum Forks ☐ With Suspension 					E1015 \$ 55
☐ With Quick Release Axles 					K0073 \$ 107
☐ Frog Leg Suspension Forks Std 1" Wide (For 3", 4" x 1" width Caster and 5" lighted Caster) (Not available with Quick Release Axles)					
Color ☐ Anodized Black ☐ Anodized Silver 					E1015 \$ 575
☐ Frog Leg Suspension Forks 1.5" Wide (For 3", 4" x 1.5" width Caster and 5" Caster) (Not available with Quick Release Axles)					
Color ☐ Anodized Black ☐ Anodized Silver 					E1015 \$ 575
Sub Total Section C 1/2					\$
Sub Total Section C 2/2					\$
Total Section C					\$

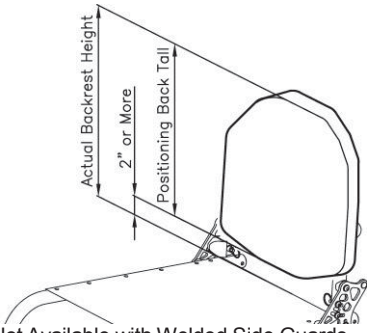
D COMPONENTS

WHEEL LOCKS								
☐ Push to Lock / Anodized Black Only		\$ STD						
☐ NEW Compact Push to Lock / Anodized Black Only ☐ w/ 2" Extension Handle (\$65up)		\$ 110						
☐ NEW Compact Scissor Lock / Anodized Black Only		\$ 110						
☐ D's Locks (*Not available with X-Core, Topolino WX2.5, Pro Wheels) (*Does not include installation fee)		\$ CFP						
PUSH HANDLES *If you are ordering push handles, You must select the backrest locking mechanism.								
☐ Integrated w/Fixed Backrest Height / Anodized Black Only		\$ STD						
☐ Bolt-on / Anodized Black Only Specify Height from floor -()		\$ 142						
☐ Bolt-on w/Quick Release / Anodized Black Only Specify Height from floor -()		\$ 190						
☐ NEW Height Adj. / Anodized Black Only Specify Maximum height from floor -()		\$ 295						
NEW ARM RESTS Anodized Black Only,								
☐ Fixed Tubular Armrest	Adj. Angle and Height Height from seat sling Set at : ☐ 8" ☐ 9" ☐ 10" ☐ 11" ☐ 12" E0973	\$ 250						
☐ Swing Away Tubular Armrest	Adj. Angle and Height Height from seat sling Set at : ☐ 8" ☐ 9" ☐ 10" ☐ 11" ☐ 12" E0973	\$ 250						
☐ Flip-Up Armrest	☐ w/Fixed Tubular Armrest Height from seat sling Set at : ☐ 8" ☐ 9" ☐ 10" ☐ 11" ☐ 12" E0973	\$ 250						
Adj. Angle ☐ w/ Swing Away Tubular Armrest	Height from seat sling Set at : ☐ 8" ☐ 9" ☐ 10" ☐ 11" ☐ 12" E0973	\$ 250						
☐ Custom Length & Height (\$150 up charge) *Specify Height & Length → 								
ANTI-TIPPERS								
☐ Anti-Tippers (pair) w/ Quick Release / Anodized Black Only		E0971 \$ 112						
SIDE GUARDS								
☐ Rigid Side Guards - Black ABS	☐ Non-Removable ☐ Removable ☐ Removable w/Hole (\$10 up charge) Size: ☐ STD ☐ Medium ☐ Cut to Wheel (\$98 up charge) Measurement Color: Black Only	\$ 142						
☐ Rigid Side Guards - Aluminum  <table border="1" data-bbox="186 1165 414 1249"> <tr> <th>Height</th> <th>Size</th> </tr> <tr> <td>7.5"</td> <td>STD</td> </tr> <tr> <td>8.25"</td> <td>Medium</td> </tr> </table>	Height	Size	7.5"	STD	8.25"	Medium	☐ Non-Removable ☐ Removable Size: ☐ Medium Only ☐ Cut to Wheel (\$98 up charge) Measurement Color ☐ Black (STD) ☐ Pol. Aluminum (STD) ☐ Powder Coated () (\$50 up charge) ☐ Gold Plate (+CFP)	\$ 197
Height	Size							
7.5"	STD							
8.25"	Medium							
								
☐ PRO-TEX Fender Guards / (Black ABS -Only) ☐ Cut to Fit Specs ☐ Do Not Install		\$ 175						
☐ Fabric Clothing Guards Black Nylon w/Velcro		\$ 78						
SEAT CUSHION -*Please see section E for other cushion choices.-								
☐ Foam Cushion w/Nylon Black Cover Thickness ☐ 1" ☐ 2" ☐ 3" ☐ 4"		\$ 78						
☐ Custom Cushion Cover (Std. Foam Cushion Only)-Select Fabric ()								
MISCELLANEOUS								
☐ Frame Protectors / Black Leather/Zipper Type		\$ 78						
☐ K0065 Spoke Protectors / Black Only ☐ 24" ☐ 25" ☐ 26" K0065		\$ 107						
☐ E0978 Seat Belt w/Auto Buckle (BP) 1.5" Black Belt Only E0978		\$ 69						
☐ Tie Down Brackets 4pcs.		\$ 210						
SPECIALTY PARTS ☐ Wheelchair Tray Size: ☐ W14"x D12" ☐ W16"x D14" ☐ W18"x D14" Color: ☐ Raw Alum. ☐ Black Alum. ☐ Polish Alum. (\$250 up charge)								
		\$ 250						
UPHOLSTERY CHOICE & COLOR *Positioning Back Available. Please see Section E								
☐ Standard Adjustable Ht Backrest Black Only		\$ STD						
☐ Adjustable Tension Backrest Black Only (*Not available with 8" Backrest Height)		\$ 135						
Custom Options	☐ Fabrics for Fun ☐ 001 ☐ 002 ☐ 003 ☐ 004 ☐ 006 ☐ 010 ☐ 011 . . . (\$95 up charge) ☐ For Std. or Adj. Style ☐ 014 ☐ 015 ☐ 016 ☐ 017 ☐ 018 ☐ 020 ☐ 021 ☐ 022 ☐ 023 ☐ 024 ☐ 025 ☐ 026 ☐ 027 ☐ 028 ☐ Backrest Leather Trim (outside panels) (\$77 up charge)							
☐ Calf Strap ☐ Black ☐ Matching Fabric (\$55 up charge) -Select Fabric() K0038		\$ 24						
☐ Matching Fabric on Padded Side Guards -Select Fabric ()		\$ 210						
Total Section D		\$						

CHAIR OPTIONS 1/2

SEAT CUSHION											
SUPRACOR® WHEELCHAIR CUSHIONS											
Stimulite® Contoured Cushion *Black Upholstery only		☐ W14"xD16" (CD1416)						E2603	\$ 525		
		☐ W16"xD16" (CD1616)						E2603	\$ 525		
		☐ W16"xD18" (CD1618)						E2603	\$ 525		
		☐ W18"xD16" (CD1816)						E2603	\$ 525		
		☐ W18"xD18" (CD1818)						E2603	\$ 525		
		☐ W18"xD20" (CD1820)						E2603	\$ 540		
		☐ W20"xD16" (CD2016)						E2603	\$ 540		
		☐ W20"xD18" (CD2018)						E2603	\$ 540		
		☐ W20"xD20" (CD2020)						E2603	\$ 555		
				☐ Extra Sensitive (\$30 up charge) ☐ Sling Style Bottom (\$45 up charge) ☐ Custom Size W(") x D(") (\$- CFP)							
Stimulite® Slimline Cushion *Black Upholstery only		☐ W14"xD16" (SM1416)						E2603	\$ 495		
		☐ W16"xD16" (SM1616)						E2603	\$ 495		
		☐ W16"xD18" (SM1618)						E2603	\$ 495		
		☐ W18"xD18" (SM1818)						E2603	\$ 495		
		☐ W18"xD20" (SM1820)						E2603	\$ 495		
		☐ W20"xD16" (SM2016)						E2603	\$ 510		
		☐ W20"xD18" (SM2018)						E2603	\$ 510		
		☐ W20"xD20" (SM2020)						E2603	\$ 510		
				☐ Extra Sensitive (\$30 up charge) ☐ Sling Style Bottom (\$45 up charge) ☐ Custom Size W(") x D(") (\$- CFP)							
		Stimulite® Corbee Pediatric *Black Upholstery only		☐ W10"xD10" (PD1010)						E2603	\$ 450
☐ W12"xD12" (PD1212)						E2603	\$ 450				
☐ W14"xD14" (PD1414)						E2603	\$ 450				
☐ Extra Sensitive (\$30 up charge) ☐ Sling Style Bottom (\$45 up charge) ☐ Custom Size W(") x D(") (\$- CFP)											
ROHO WHEELCHAIR CUSHIONS											
Choose Size	☐ 4" HIGH PROFILE® QUADTRO SELECT® CUSHION *Black Upholstery							E2622-2623	\$ 525		
	☐ W12"xD12"	☐ W12"xD13"	☐ W13"xD13"	☐ W13"xD15"	☐ W13"xD16"	☐ W15"xD13"	☐ W15"xD15"				
	☐ W15"xD16"	☐ W15"xD18"	☐ W16"xD15"	☐ W16"xD16"	☐ W16"xD18"	☐ W16"xD20"	☐ W18"xD15"				
	☐ W18"xD16"	☐ W18"xD18"	☐ W18"xD20"	☐ W18"xD22"	☐ W20"xD16"	☐ W20"xD18"	☐ W20"xD20"				
	☐ W20"xD22"	☐ W22"xD18"	☐ W22"xD20"	☐ W22"xD22"	☐ Custom Size W(") x D(") (\$+85)						
Choose Size	☐ 3" MID PROFILE® QUADTRO SELECT® CUSHION *Black Upholstery only							E2622-2623	\$ 525		
	☐ W15"xD15"	☐ W15"xD16"	☐ W16"xD16"	☐ W16"xD18"	☐ W18"xD16"	☐ W18"xD18"	☐ Custom Size W(") x D(") (\$+85)				
Choose Size	☐ 2" LOW PROFILE® QUADTRO SELECT® CUSHION *Black Upholstery only							E2622-2623	\$ 525		
	☐ W12"xD12"	☐ W12"xD13"	☐ W13"xD13"	☐ W13"xD15"	☐ W13"xD16"	☐ W15"xD13"	☐ W15"xD15"				
	☐ W15"xD16"	☐ W15"xD18"	☐ W16"xD15"	☐ W16"xD16"	☐ W16"xD18"	☐ W16"xD20"	☐ W18"xD15"				
	☐ W18"xD16"	☐ W18"xD18"	☐ W18"xD20"	☐ W20"xD16"	☐ W20"xD18"	☐ W20"xD20"	☐ Custom Size W(") x D(") (\$+85)				
Choose Size	☐ 2" NEXUS SPIRIT® CUSHION *Black Upholstery only							E2622	\$ 525		
	☐ W14"xD14"	☐ W14"xD16"	☐ W15"xD15"	☐ W16"xD16"	☐ W16"xD18"	☐ W17"xD17"	☐ W18"xD16"				
☐ * Custom Size Not Available											
STAR CUSHIONS											
Galaxy Cushion *Black Upholstery only		☐ W14"xD14" (GXY 1414-2) ☐ W16"xD14" (GXY 1614-2) ☐ W16"xD18" (GXY 1618-2) ☐ W18"xD18" (GXY 1818-2) ☐ W20"xD16" (GXY 2016-2) ☐ W20"xD20" (GXY 2020-2) ☐ W14"xD16" (GXY 1416-2) ☐ W16"xD16" (GXY 1616-2) ☐ W18"xD16" (GXY 1816-2) ☐ W18"xD20" (GXY 1820-2) ☐ W20"xD18" (GXY 2018-2)					Choose Zoning Options ☐ Standard - Dual Valve ☐ Dual Valve - Front to Back ☐ Four Valve ☐ Simple - Single Valve ☐ Four Valve - Enhanced ☐ Four Valve - Front to Back			E2624	\$ 498
THE CURVE CUSHION *Black Upholstery only		☐ W10"xD10"	☐ W10"xD11"	☐ W10"xD11"	☐ W10"xD12"	☐ W11"xD12"	☐ W12"xD12"	\$ 133			
		☐ W12"xD13"	☐ W12"xD14"	☐ W13"xD11"	☐ W13"xD12"	☐ W13"xD13"	☐ W13"xD14"				
		☐ W14"xD14"	☐ W14"xD15"	☐ W14"xD16"	☐ W15"xD14"	☐ W15"xD15"	☐ W15"xD16"				
		☐ W16"xD14"	☐ W16"xD15"	☐ W16"xD16"	☐ W16"xD17"	☐ W16"xD18"	☐ W17"xD16"				
		☐ W17"xD17"	☐ W17"xD18"	☐ W18"xD16"	☐ W18"xD17"	☐ W18"xD18"					
NOTE											

E CHAIR OPTIONS 2/2

POSITIONING BACK			
*Note : Actual height of the backrest will increase 2" or more.			
ADI CARBON BACKREST (Support Contour-ACTIVE 2.5") *Not Available with Welded Side Guards			
Mounting Hardware: ☐ 2 Point Pro ☐ Quick Release			
Height	☐ LOW (10" Tall) *Backrest Height requires 12" or more	*Available 13"-19" width	Width () E2613 \$ 745
	☐ Medium (13" Tall) *Backrest Height requires 15" or more	*Available 13"-19" width	Width () E2613 \$ 815
	☐ Standard (16" Tall) *Backrest Height requires 18" or more	*Available 13"-19" width	Width () E2613 \$ 880
	☐ Custom Size (" Tall) x (" width)		E2613 \$ CFP
ADI ALUMINUM BACKREST (Support Contour-ACTIVE 2.5") *Not Available with Welded Side Guards			
Mounting Hardware: ☐ 2 Point Pro ☐ Quick Release			
Height	☐ LOW (10" Tall) *Backrest Height requires 12" or more	*Available 10"-20" width	Width () E2613 \$ 610
	☐ Medium (13" Tall) *Backrest Height requires 15" or more	*Available 12"-20" width	Width () E2613 \$ 625
	☐ Standard (16" Tall) *Backrest Height requires 18" or more	*Available 14"-20" width	Width () E2613 \$ 640
	☐ Tall (18" Tall) *Backrest Height requires 20" or more	*Available 16"-20" width	Width () E2613 \$ 655
	☐ X-Tall (20" Tall) *Backrest Height requires 22" or more	*Available 16"-20" width	Width () E2613 \$ 670
	☐ Custom Size (" Tall) x (" width)		E2613 \$ CFP
ADI ALUMINUM BACKREST (Support Contour-DEEP 4") *Not Available with Welded Side Guards *Not Available with Lateral Supports			
Mounting Hardware: ☐ 2 Point Pro ☐ Quick Release			
Height	☐ LOW (10" Tall) *Backrest Height requires 12" or more	*Available 10"-20" width	Width () E2615 \$ 710
	☐ Medium (13" Tall) *Backrest Height requires 15" or more	*Available 12"-20" width	Width () E2615 \$ 725
	☐ Standard (16" Tall) *Backrest Height requires 18" or more	*Available 14"-20" width	Width () E2615 \$ 740
	☐ Tall (18" Tall) *Backrest Height requires 20" or more	*Available 16"-20" width	Width () E2615 \$ 755
	☐ X-Tall (20" Tall) *Backrest Height requires 22" or more	*Available 16"-20" width	Width () E2615 \$ 770
	☐ Custom Size (" Tall) x (" width)		E2615 \$ CFP
ADI LATERAL SUPPORTS (Lateral Supports are compatible with only active 2.5" ADI positioning back.)			
☐ Stationary Pad	Pad Size ☐ Small (3" H x 3" D) ☐ Medium(3" H x 4" D) ☐ Large(4.5" H x 4.5" D) ☐ Large(4" H x 6" D)		
	☐ Pair ☐ Left side only or ☐ Right side only		
☐ Swing-Out	Pad Size ☐ Small (3" H x 3" D) ☐ Medium(3" H x 4" D) ☐ Large(4.5" H x 4.5" D) ☐ Large(4" H x 6" D)		
	☐ Pair ☐ Left side only or ☐ Right side only		
OTHER			
☐ ADI Transfer Board ☐ 21" (ADV21) ☐ 24" (ADV24) ☐ 29"			
☐ ADI Push and Transfer Gloves Leather (one size fits)			
☐ ADI Carbon Solid Seat (Wide: 18" Max, Depth: 20") ※Colours Side Guards Not Available with ADI Carbon Solid Seat			
With / ☐ ADI Folding Carbon Side Guard			
With / ☐ ADI Fixed Carbon Side Guard			
VARILITE® Seating Systems Available		Please visit for all options. www.varilite.com	
NOTE			
		Sub Total Section E 1/2 \$	
		Sub Total Section E 2/2 \$	
		Total Section E \$	

F SPECIAL PARTS COLORS (Please allow an additional 1-2 weeks delivery)
PARTS COLORS

Powder-Coated	Gloss Black	Twilight Flake	Yellow	Magenta	Fire Engine Red	White
	Candy Blue	Candy Purple	Candy Green	Candy Red	Candy Raspberry	Cobalt Blue
	Matte Gray	Matte Purple	Matte Blue	Matte Green	Matte Black	Matte Red
	Neon Red	Neon Green	Orange Sparkle	Mellon Yellow	Gun Smoke Silver	Pigeon Blue
		Granny Sparkle Green	Sparkle Sky Blue	Vermilion	Reed Green	Blue Lilac

Special Colors Polished Aluminum
24 karat Gold Plated

Please Choose the Color for each Part

*Parts not selected will be black color. (Suspension Spring: Silver Color)

*Minimum 3 parts selected per order.

*PC=Powder-Coated

Saddle Stops

☐ PC ()	\$40
☐ Polished	\$115
☐ Gold Plated	\$CFP

Backrest

☐ PC ()	\$125
☐ Polished	\$225
☐ Gold Plated	\$CFP

Backrest Brackets

☐ PC ()	\$90
☐ Polished	\$115
☐ Gold Plated	\$CFP

Foot Plate

☐ PC ()	\$45
☐ Polished	\$90
☐ Gold Plated	\$CFP

Footrest Rings

☐ PC ()	\$20
☐ Polished	\$40
☐ Gold Plated	\$CFP

Foot Plate Clamps

☐ PC ()	\$50
☐ Polished	\$75
☐ Gold Plated	\$CFP

Footrest Tubes

☐ PC ()	\$65
☐ Polished	\$85
☐ Gold Plated	\$CFP

Caps

☐ PC ()	\$20
☐ Polished	N/A
☐ Gold Plated	N/A

Caster Housing

☐ PC ()	\$150
☐ Polished	\$280
☐ Gold Plated	N/A

Front Springs

☐ PC ()	\$20
☐ Polished	N/A
☐ Gold Plated	N/A

Forks & Knurl Nuts

☐ PC ()	\$65
☐ Polished	\$90
☐ Gold Plated	\$CFP

Axle Plates

☐ PC ()	\$65
☐ Polished	\$85
☐ Gold Plated	\$CFP

Camber Plates

☐ PC ()	\$50
☐ Polished	\$75
☐ Gold Plated	\$CFP

Axle Bases

☐ PC ()	\$50
☐ Polished	\$75
☐ Gold Plated	\$CFP

Axle Clamps

☐ PC ()	\$50
☐ Polished	\$75
☐ Gold Plated	\$CFP

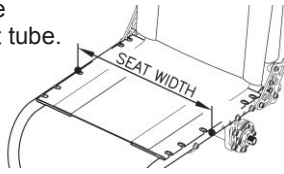
Axle Tube

☐ PC ()	\$50
☐ Polished	\$75
☐ Gold Plated	\$CFP

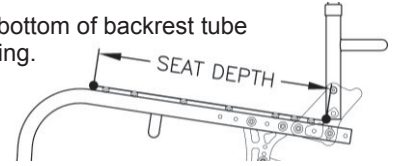
Total Section F \$

1. SEAT WIDTH

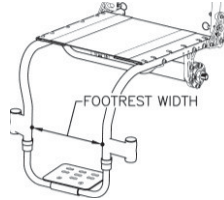
Measure from outside of seat tube to the outside of the opposite seat tube.

**2. SEAT DEPTH**

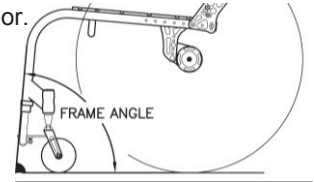
Measure from at the front bottom of backrest tube to the front edge of seat sling.

**3. FOOTREST WIDTH**

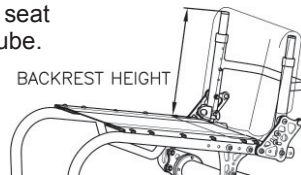
Measure from inside of front frame tube to inside of opposite front frame tube.

**4. FRAME ANGLE**

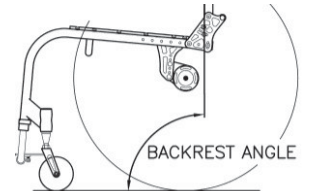
Measure angle frame front to floor.

**5. BACKREST HEIGHT**

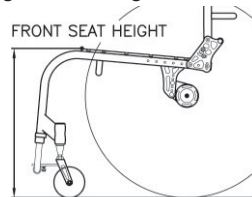
Measure from top of back post to seat sling at front bottom of backrest tube.

**6. BACKREST ANGLE**

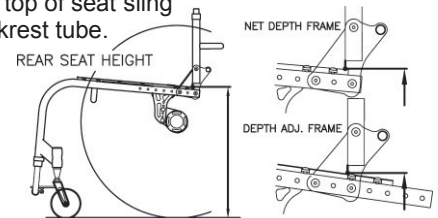
Measure angle from front of back post to floor.

**7. FRONT SEAT HEIGHT**

Measure from floor to top of seat sling at front edge.

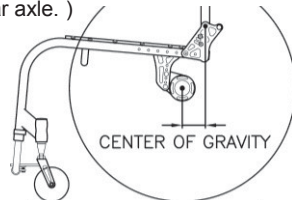
**8. REAR SEAT HEIGHT**

Measure from floor to top of seat sling at front bottom of backrest tube.

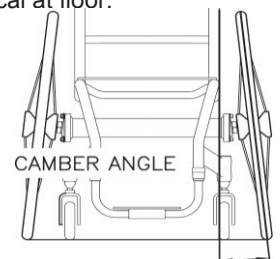
**9. CENTER OF GRAVITY**

Measure from center of back tube at the folding pivot to the center of the rear axle.

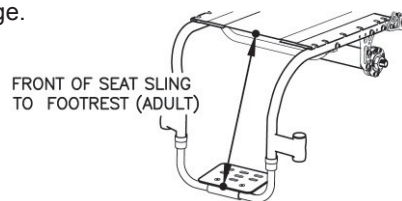
(If welding backrest is selected, Measure from center of back tube at the seat tube to the center of the rear axle.)

**10. CAMBER ANGLE**

Measure angle from wheel to vertical at floor.

**11. FRONT OF SEAT SLING TO FOOTREST****ADULT CHAIR**

Measure from top of seat sling at front edge to top of foot plate at front edge.

**PEDIATRIC CHAIR**

Measure vertical from top of seat sling at front edge to top of foot plate.

**12. INSIDE WHEEL CLEARANCE**

Measure from outside of seat back tube to inside of rear tire.

(If side guards are selected, measure from outside of side guards to inside of rear tire.)

