

SUMMIT SL3500D


OUTDOOR CONVENIENCE



The SL3500D outdoor stairlift is engineered with a durable extruded aluminum track and marine-grade vinyl seating to handle rain, wind, snow, and temperatures from 0°F to 125°F. Includes header, internal chassis, and a large weatherproof cover for protection from outside elements.

FEATURES AND BENEFITS

- Designed for long-term outdoor use
- Header cover, internal chassis cover, and a large weatherproof cover included
- Rocker controller in armrest operates with constant-pressure functionality to stop the lift immediately if the control is released
- Extruded aluminum rail mounts to the stairs (not the wall) on either side of staircase
- Operates on standard 115V AC household current; no batteries required
- Can be installed at 350 lb capacity up to 15'9" and 250 lb up to 25'
- Call/send controls at top and bottom of stairs

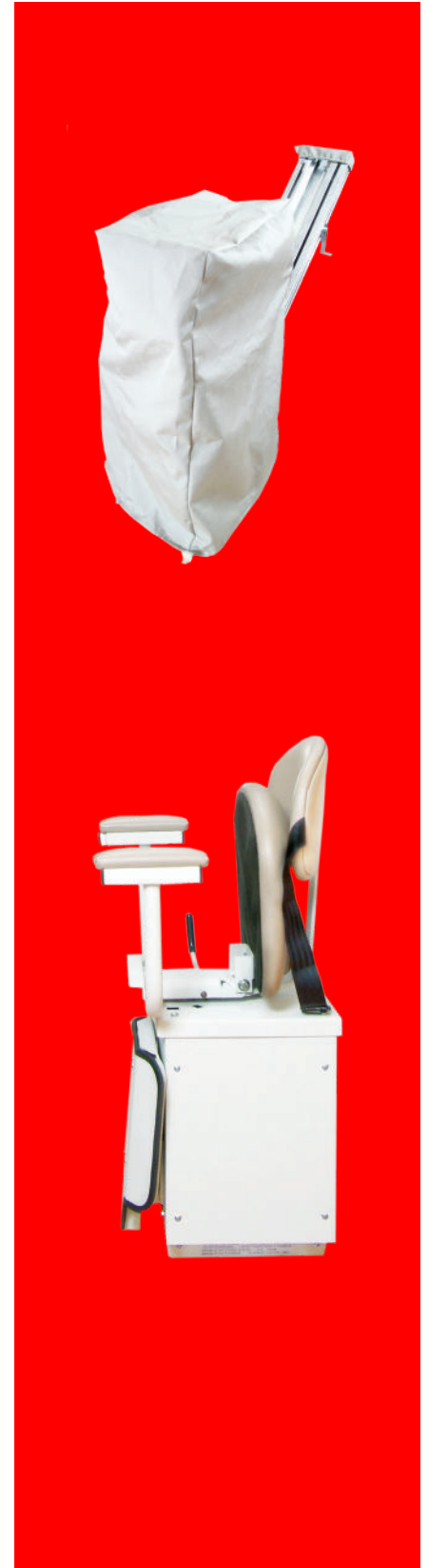

Constant pressure
rocker control

Seat belt



Footrest obstruction switch



My wife had trouble walking up to our front door from the street. Having the SL3500D allows us to enjoy our forever home for years to come by giving her access to the home year around.

SL3500D PRODUCT SPECIFICATIONS

WEIGHT CAPACITY - 350 lb up to 15' 9" track
- 250 lb for 15'10" - 25' track

TRACK LENGTH 16' to 25'

SEAT SIZE 17.5" x 14"

FOLDED WIDTH 14" from wall

**SAFETY
FEATURES**

- Footrest obstruction switch
- Swivel seat cut-off switch
- Final limits
- Slack cable brake with switch
- Constant pressure controls
- Seat belt

INCLINE LIMITS 35° to 45°

SPEED 20 FPM

FOOTREST SIZE 9.5" x 10"

DRIVE Cable

WARRANTY 3-year parts

- OPTIONS**
- Key lock
 - Extra rail